

Power assist device

SmartDrive MX2+



permobil

On-demand power assist

Pushing a manual wheelchair can put stress on your shoulders and lead to serious future issues. SmartDrive's seamless power assist has been proven to reduce up to 81%* of your pushes and minimise repetitive stress, allowing you to focus on the things that matter.

*based on data collected from SmartDrive's SmartEvaluation App



Multiple ergonomic handles to quickly attach to almost any manual wheelchair



Easily removable, making it great for transportation



Various control options to power your push, your way

OmniWheel's patented design engages and stays engaged in any position, maintains optimal traction, and climbs up to 10 degree slope (17% inclines)

Lightweight power assist weighing in at 6kg

SwitchControl

SwitchControl is a simple way to control your SmartDrive with the push of a button. You can use SwitchControl in momentary or latched mode and optionally, attach an external switch with a custom mono jack option.

Momentary

Your SmartDrive will engage for as long as you hold either Switch.

Latched

Press and hold either Switch to set your cruise speed. When you want to stop, simply tap the Switch and hold on to your wheels.

SwitchControl with Mono Jack Option

Use any external mono jack switch to operate your SwitchControl; no programming required.



**Buddy Button with
extended mount and
SwitchControl with
mono jack**




Premium wearable controllers

Complement your wired controller with a premium wearable device!

You can use your own Apple Watch or the Samsung Galaxy Watch to control your SmartDrive.

- Fast and simple setup
- Increased screen visibility
- Improved battery performance
- Free app from App Store or Google Play Store
- Refer to our website for full list of compatible wearables

* Apple Watch, App Store, Samsung Galaxy Watch, and Google Play Store are trademarks or registered trademarks of their respective owners.

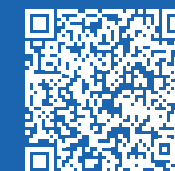

Scan to learn more about
SmartDrive MX2+
Premium Wearables

Power your push, your way

Discover how our variety of hardware and software solutions can be combined to better support your individualised mobility goals.



Specifications

Dimensions	Length: 389mm (15.3") Width: 141mm (5.7") Height: 242mm (9.5")
Product weight	6kg (13.2lb) Drive Unit
User weight	14kg to 150 kg (30lb to 331 lb)
Driving range*	Up to 19.8 km (12.3 mi)
Motor power	250W Brushless DC
Maximum speed	8.8km/h (5.5mph)
Operating temp	-25°C to 50°C (-13°F to 122°F)
Chair type compatibility	Rigid, Folding, Tilt-In Space, One-Arm Drive and others Chair widths: 8"(no anti-tips) -22"
Drive wheel diameter	501mm, 540mm, 559mm, 590mm (22", 24", 25", 26")
Warranty	2-Year Limited Warranty
Bluetooth LE	FCC ID: QOQBT113 IC ID: 5123A-BGTBLE113

**Performance measurements will vary based on user and chair characteristics, driving and battery conditions.*

The SmartDrive has been tested and conforms to all applicable requirements of ANSI/RESNA Standards for Wheelchairs - Volume 1 and Volume 2 and EN 12184.

Note: A wheelchair's specifications provided by its manufacturer may be slightly affected by the addition of the SmartDrive. Please refer to your owner's manual for more detailed operation and care instructions.