

Power wheelchairs

# F-Series



permobil

# Comfort is king

Corpus®, our most ergonomic and capable power wheelchair seating system, was designed to configure to your body. Biomechanically optimised power positioning enables functional movements while supporting your optimal position throughout the day.

Dual density backrest foam promotes greater surface area contact while the contoured Ergo seat pelvic well insert can help pressure redistribution and reduce peak pressures. STRETCH-AIR® covers offer a breathable, moisture-resistant fabric material to enhance comfort and help protect underlying foam.

STRETCH-AIR  
breathable covers



180° of legrest elevation

180° of power recline  
50° of posterior tilt

# BodiLink on Corpus

Specially designed for the Corpus seating system, BodiLink accessories feature best-in-class hardware. Offering true anatomical support with the highest level of adjustability and comfort to meet your individual positioning needs.

Head Support—  
Easy to adjust piano hinge and double-taper joints for maximum strength

Lateral Pelvic/Thigh Support— Easily position along the UniTrack rail with 360° pad adjustment and user-friendly quick release lever



## Inspired design

The F-Series builds on over 50 years of experience to create our most capable family of front-wheel drive power wheelchairs ever. We've combined a low center of gravity, optimised weight distribution and powerful high efficiency motors to offer enhanced base stability, obstacle climbing and chair responsiveness.

**ComfortRide Pro™**—  
Oil dampened shocks for  
greater vibration isolation

**CasterGlide™**—  
Worry free, adjustment  
free, anti-flutter device

**ComfortRide**

F3 Corpus

**ComfortRide Pro**

F5 Corpus & F5 Corpus VS

## Superior suspension

Conquer your world with ComfortRide™ and ComfortRide Pro™ fully independent, adjustable suspension set to meet your needs. Enjoy superior vibration reduction, ride comfort and drive performance for the widest variety of driving scenarios. The front-wheel driving experience is further enhanced by our patented ESP (Enhanced Steering Performance) gyroscope tracking technology, which provides exceptional tracking at high speeds, on side slopes or with alternative drive systems.



3.5" rearward shift during elevate reduces footprint to get you closer to tasks

14"

3 points of contact for added stability

## Reach for more

ActiveHeight™ and ActiveReach™ combine seamlessly to give you the ultimate advantage in accessing your environment. Seat elevation with integrated forward tilt creates the most innovative, functional and stable seating

solution on the market today. With up to 14" of vertical seat adjustment and 45° of ActiveReach forward tilt capability, the F-Series empowers you to reach for more and go further faster with new elevated drive speeds of up to 5.1 km/h.

	F3	F5	F5 VS
↑ <b>ActiveHeight</b>	12"	14"	14"
↶ <b>ActiveReach</b>	30°	45°	Standing

*10° of no charge ActiveReach with ActiveHeight*



## Small footprint, big impact

The F3 Corpus offers the smallest footprint of the F-Series front-wheel drive power wheelchairs. Its compact size and agility allows users to easily manoeuvre in and out of tight spaces. The optimised battery weight distribution, lower center of mass and fully adjustable ComfortRide suspension improve stability while traversing inclines or declines. Experience the enhanced functionality of the F3 by combining outstanding chassis performance with Permobil's industry-leading Corpus Seating System offering 12" of ActiveHeight and up to 30° of ActiveReach forward tilt technology.

F3 Corpus



36" without support wheels

## Strength and stability

Power meets performance on the F5 Corpus with a 12 km/h top end speed and upgraded ComfortRide Pro suspension. All-wheel independent, fully adjustable, oil dampened shocks provide the ultimate in ride comfort while helping to reduce vibrations that can transfer from terrain up to user. Experience enhanced vibration reduction driving in all seat positions including an unprecedented 14" ActiveHeight seat elevation. Climb and reach higher with increased obstacle negotiation, exceptional motor torque and enhanced stability allowing up to 45° of ActiveReach forward tilt functionality. Standard performance options include aggressive tread tyres and castor mud guard protection.

F5 Corpus



Aqua

## Stand without compromise

We combined all of the performance features of the F5 Corpus with superior power standing functionality. Fully programmable standing sequences allow user position to be optimised using supine, semi-reclined or sit-to-stand techniques. Further customise standing positions with adjustable chest support, mouldable knee blocks and power articulating footplates. From face-to-face interactions to improved health outcomes, the F5 Corpus VS empowers you to stand with confidence.

Adjustable chest support

F5 Corpus VS



Gold

High power LED headlights are twice as bright, making it easier to see and be seen



Metallic colours

Make a statement

With a sleek design, expanded colour/hydrographic options and powerful LED lighting, you should be prepared to turn a few heads in your new F-Series.

Please note: Actual product colours may vary from colours shown



Gold




Platinum



Bronze

Hydrographic options

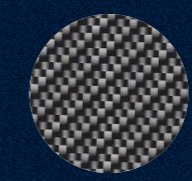


Carbon Fiber



Mossy Oak® Breakup® Country™



Water Lily

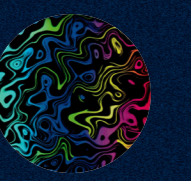

Liquid Fusion\*

\*Liquid Fusion is currently only available on the Corpus VS models

## Specifications



	F3 Corpus	F5 Corpus	F5 Corpus VS
<b>Maximum user weight</b>	150 kg	150 kg	136 kg
<b>Maximum speed</b>	10 km/h	12 km/h	12 km/h
<b>Range</b>	32.2 km <sup>*</sup>	32.2 km <sup>*</sup>	32.2 km <sup>*</sup>
<b>Base width</b>	24"	25.5"	25.5"
<b>Base length</b>	40" with support wheels 36" without support wheels	42" with support wheels 38" without support wheels	43"
<b>Base turning radius</b>	26.5"	30"	30"
<b>Ground clearance</b>	3"	3"	3"
<b>Recommended battery type</b>	Group 24	Group 24	Group 24
<b>Weight including batteries</b>	185 kg	187 kg	196 kg
<b>Drive electronics</b>	R-net 120A	R-net 120A	R-net Power Platform 120A
<b>Seat-to-floor height</b>	17.5", 18.5", 19.5" <sup>†</sup>	17.5", 18.5", 19.5" <sup>†</sup>	17.5", 18.5", 19.5" <sup>†</sup>
<b>Seat elevation</b>	12"	14"	14"
<b>Drive wheel</b>	14"	14"	14"
<b>Castors</b>	8"	8"	8"
<b>Tilt options</b>	Posterior: 0° – 50° ActiveReach: 5°, 10°, 20°, 30° (up to 136 kg)	Posterior: 0° – 50° ActiveReach: 5°, 10°, 20°, 45° (up to 136 kg)	Posterior: 0° – 50° Standing
<b>Legrest angle- power</b>	85° – 180°	85° – 180°	85° – 180°
<b>Backrest angle- power recline</b>	85° – 180°	85° – 180°	85° – 180°
<b>Backrest height</b>	20", 23" – 28" (1" incr.)	20", 23" – 28" (1" incr.)	23" – 28" (1" incr.)
<b>Seat depth</b>	14" – 22" (1" incr.)	14" – 22" (1" incr.)	14" – 22" (1" incr.)
<b>Seat width</b>	17" – 23" (2" incr.)	17" – 23" (2" incr.)	17" – 23" (2" incr.)

<sup>\*</sup> Actual driving range will vary based on driving and battery conditions. | <sup>†</sup> Measured without cushion. | Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and includes base, batteries, power seat functions and seating system without accessories.



**Note:** MyPermobil connectivity is an available option on all new F-Series Corpus & Corpus VS models. Power Platform and +Connect form the core digital infrastructure of this power wheelchair, combining advanced assistive technology with smart connectivity. Connected from the factory\*, this chair is online and ready for proactive service, remote insights, and efficient support from day one. The connected chair reflects Permobil's commitment to reliability, transparency, and confidence over the life of the chair.

\*From 1 July, 2026, for chairs delivered in New Zealand.

App Store is a trademark of Apple Inc. Google Play is a trademark of Google LLC.

