

Power wheelchairs

# Power standing



permobil

# Corpus VS

Fully programmable standing sequences allow user position to be optimised using supine, semi-reclined, or sit-to-stand techniques. Further customise standing positions with adjustable chest support, multiple knee support options, and power articulating footplates. From face-to-face interactions to improved health outcomes, the Corpus VS seating system empowers you to stand with confidence.




## Designed for independence

New swing-away knee supports are 75% lighter\* and feature improved comfort and flexible on-chair storage.

*\*As compared to original single-post design*



Release mechanism designed for a wide variety of hand function

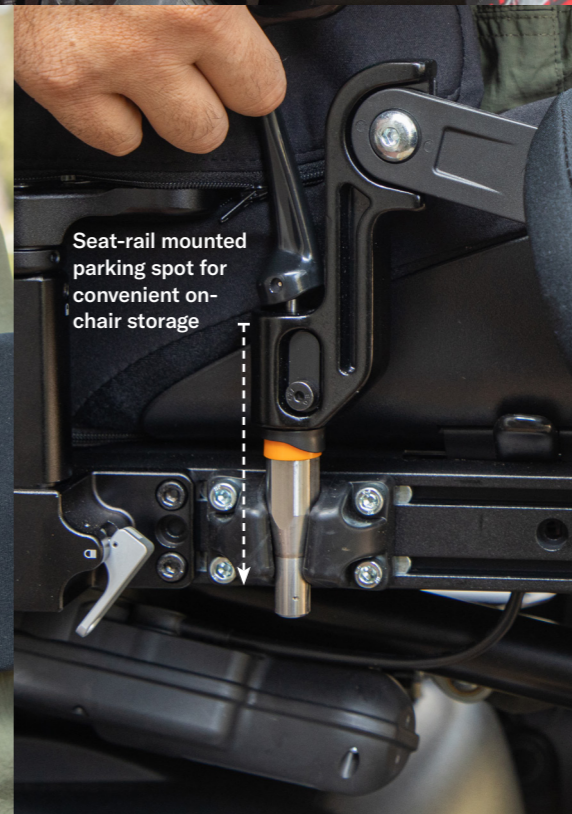

Seat-rail mounted parking spot for convenient on-chair storage



Integrated BOA® Fit System for a precision fit

*BOA is a trademark of BOA.*

- Option available on both F5 Corpus VS & M Corpus VS
- Option available only on specific chair


## Made to move

The M Corpus VS brings proven, customisable power standing to the manoeuvrability of a mid-wheel drive base. DualLink suspension balances stability, traction, and ride comfort and features fully adjustable, oil dampened shocks while excellent low-end torque ensures a smooth stand and drive experience.



## Stand without compromise

The F5 Corpus VS combines all the performance features of the F5 Corpus front-wheel drive base with superior power standing functionality. With ComfortRide Pro™ suspension, experience all-wheel independent, fully adjustable, oil dampened shocks for the ultimate in ride comfort while optimising traction and performance.

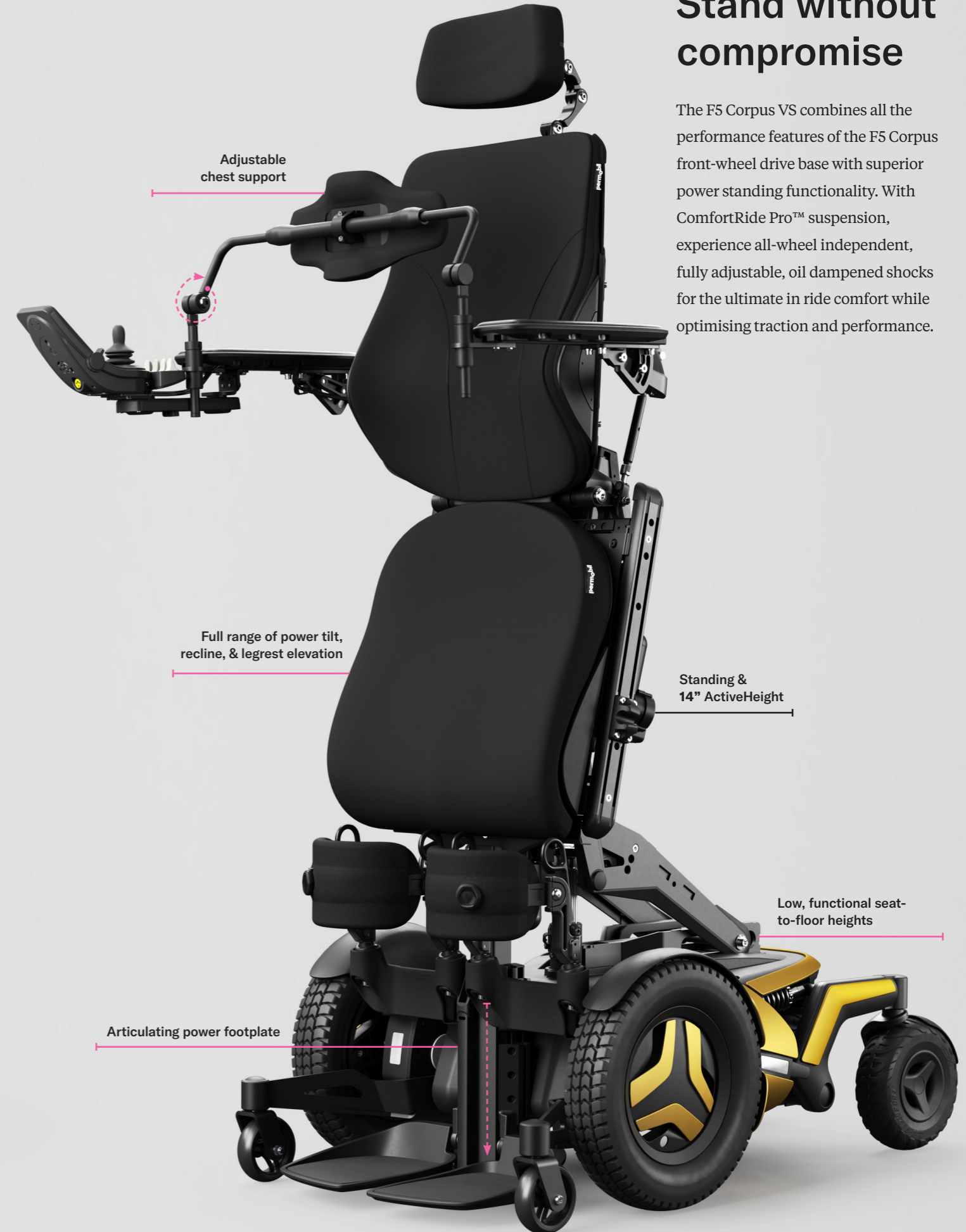

No additional hardware required

## Programming simplified

Programming changes are faster and more intuitive with QuickConfig, the new wireless programming app for Permobil wheelchairs equipped with Power Platform. Designed with clinicians and equipment providers in mind, in-app visuals ease customisation of standing sequences, memory positions, and drive performance. The result? Wheelchair setup has never been easier.



Power Platform and +Connect form the core digital infrastructure of this power wheelchair, combining advanced assistive technology with smart connectivity. Connected from the factory\*, this chair is online and ready for proactive service, remote insights, and efficient support from day one. The connected chair reflects Permobil's commitment to reliability, transparency, and confidence over the life of the chair.

\*From 1 July, 2026, for chairs delivered in New Zealand.



## Specifications

|                                       | M Corpus VS                      | F5 Corpus VS                     |
|---------------------------------------|----------------------------------|----------------------------------|
| <b>Maximum user weight</b>            | 120 kg                           | 136 kg                           |
| <b>Maximum speed</b>                  | 10 km/h                          | 12 km/h                          |
| <b>Range</b>                          | 32.2 km*                         | 32.2 km*                         |
| <b>Base width</b>                     | 25"                              | 25.5"                            |
| <b>Base length</b>                    | 36.5"                            | 43"                              |
| <b>Base turning radius</b>            | 22"                              | 30"                              |
| <b>Ground clearance</b>               | 3"                               | 3"                               |
| <b>Recommended battery type</b>       | Group 24                         | Group 24                         |
| <b>Weight including batteries</b>     | 193 kg                           | 196 kg                           |
| <b>Drive electronics</b>              | R-net Power Platform 120A        | R-net Power Platform 120A        |
| <b>Seat-to-floor height</b>           | 17.5", 18.5", 19.5" <sup>†</sup> | 17.5", 18.5", 19.5" <sup>†</sup> |
| <b>Seat elevation</b>                 | 12"                              | 14"                              |
| <b>Drive wheel</b>                    | 14"                              | 14"                              |
| <b>Castors</b>                        | 7"                               | 8"                               |
| <b>Tilt options</b>                   | Posterior: 0° – 50°              | Posterior: 0° – 50°              |
| <b>Power standing</b>                 | Up to 70°                        | Up to 80°                        |
| <b>Legrest angle - power</b>          | 85° – 180°                       | 85° – 180°                       |
| <b>Backrest angle - power recline</b> | 85° – 180°                       | 85° – 180°                       |
| <b>Backrest height</b>                | 23" – 28" (1" incr.)             | 23" – 28" (1" incr.)             |
| <b>Seat depth</b>                     | 14" – 22" (1" incr.)             | 14" – 22" (1" incr.)             |
| <b>Seat width</b>                     | 17" – 23" (2" incr.)             | 17" – 23" (2" incr.)             |

\* Actual driving range will vary based on driving and battery conditions. | <sup>†</sup> Measured without cushion. | Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions. Chair weights are approximate and include base, batteries, power seat functions, and seating system without accessories.

App Store is a trademark of Apple Inc. Google Play is a trademark of Google LLC.



Permobil New Zealand | 0800 115 222 | sales.nz@permobil.com | www.permobil.com/en-nz

BR-NZ-PS002.0