

Seating & positioning

Product catalog



permobil

Table of contents

Back supports

Cover technology	3
Hardware technology	4 - 5
Back support technology.....	6
General use back supports	7
Positioning back supports	8 - 14
Positioning back supports with laterals	15 - 16

Cushions

Cushion features	19
General use	20 - 21
Skin protection	22
Positioning	22
Skin protection and positioning.....	23
Adjustable skin protection.....	23 - 27
Adjustable skin protection and positioning	27 - 31
Custom	32 - 33
Non coded cushions	33 - 34
ROHO® cushion sizing guide	34

Supports & accessories

BodiLink®	37 - 43
Lower extremity	44 - 48
Upper extremity	49 - 52
Solid seat support	52 - 53
Replacement covers.....	53 - 54
Additional accessories.....	55 - 58


References

A side-view photograph of a person sitting in a wheelchair. The person is wearing a red short-sleeved shirt and white pants. They are using a black wheelchair with a complex back support system. The backrest is reclined, and there is a large, black, mesh-like cushion or support structure behind the person's head and neck. The wheelchair has large black wheels and a smaller front wheel. The background is a blurred outdoor setting with palm trees and a building.

Back supports

Permobil back support systems represent the leading edge in innovative, customizable wheelchair back supports. Our full line of products ensure the ideal fit for every user. From general use to a complete custom system, our back supports offer seating solutions as unique as the individual using them.

Cover technology

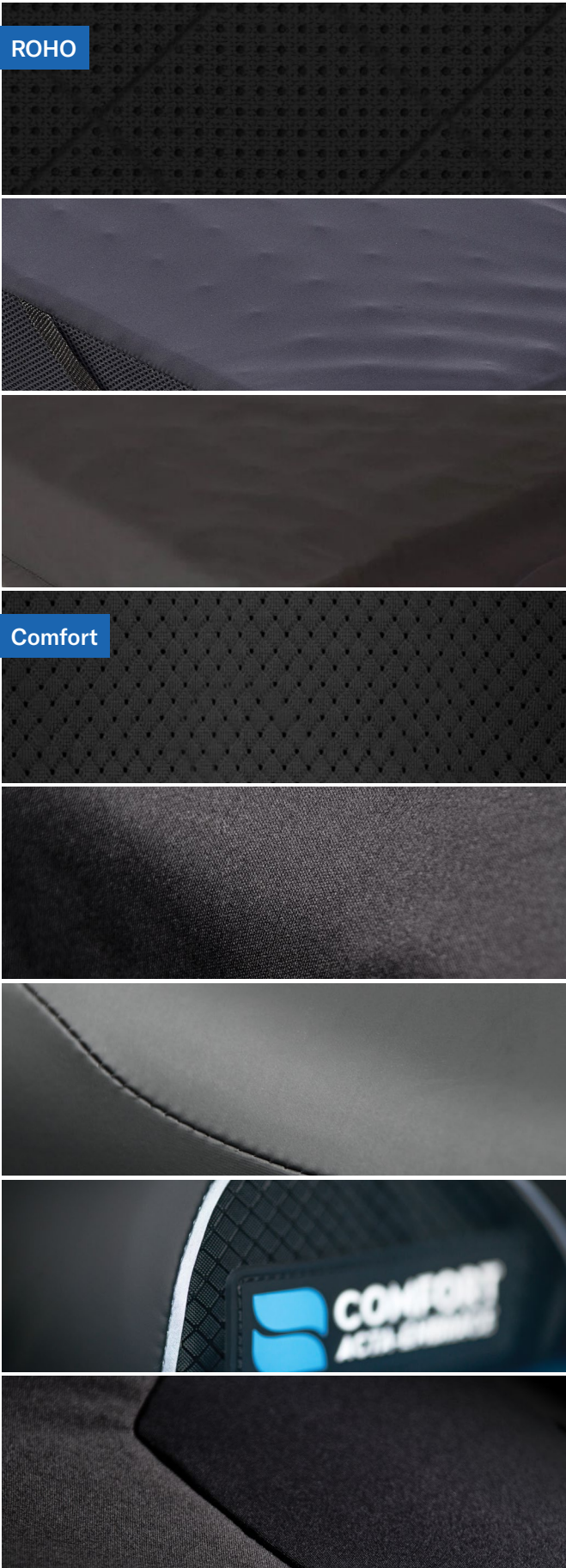

Spacer fabric

Made from durable, vapor-permeable polyester for superior breathability.

Standard cover

Multidirectional stretch fabric designed for comfort and optimal immersion.

Heavy duty cover

Multidirectional stretch, fluid-resistant fabric for protection against incontinence.

Diamond Spacer

3mm thick, 100% polyester diamond spacer fabric is designed to be breathable and allow for optimal immersion.

Stretch-Air

Our Stretch-Air is designed to be breathable with multi-directional stretch for optimal immersion.

Comfort-Tek

Combines strength and softness with multi-directional stretch for a custom fit, while providing an easy-to-clean, infection-control and incontinence barrier.

Reflective piping

Comfort back supports and cushions feature reflective piping that returns light to its source, enhancing visibility and promoting safety.

GlideWear™

An optional upgrade for Comfort cushions and Comfort Unite Back Supports and BodiLink accessories, GlideWear technology reduces pressure injury risk while enhancing comfort.

Note: All cover options shown are latex free.

Hardware technology

Compass® 4 hardware


Compass 4 hardware is strong, light weight, and independently adjustable in depth and angle of recline.

Features



Fit

Cane clamp
Fits expanded size range of 3/4", 7/8", 1", and 1 1/8" diameter tubing, and the clamp only requires 1 1/2" of space on each cane for mounting.
Industry leading rotational clamp strength
Compass 4 = 250 lbs

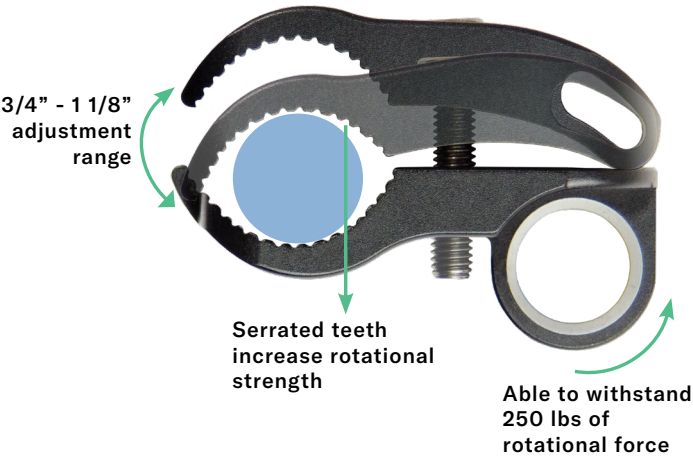


Adjustability

Range of depth adjustability


- Standard configuration (as shipped)
- Flipped cane clamp
- Flipped cane clamp & angle/depth bracket
- Flipped 90° bracket & cane clamp
- Flipped 90° bracket
- Cane tube center line

Each configuration has 2" of depth adjustment
Showing depth range from center of cane clamp to the support surface (1.5" thick)

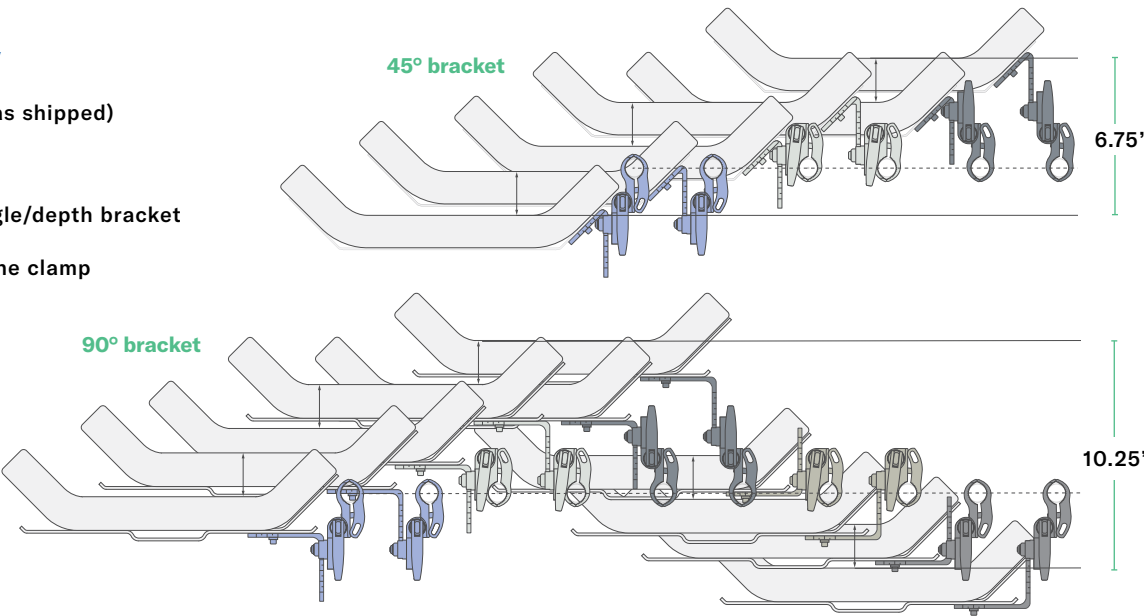

Adjustability

Width adjustability

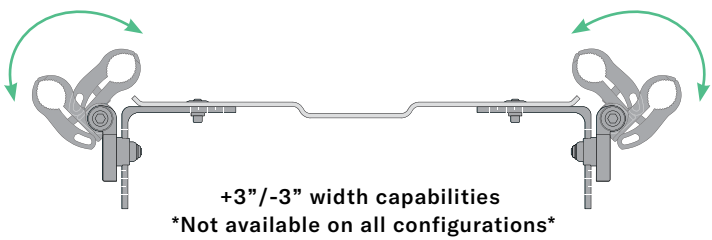


Angle adjustability



Depth adjustability



Isolated depth & recline adjustability

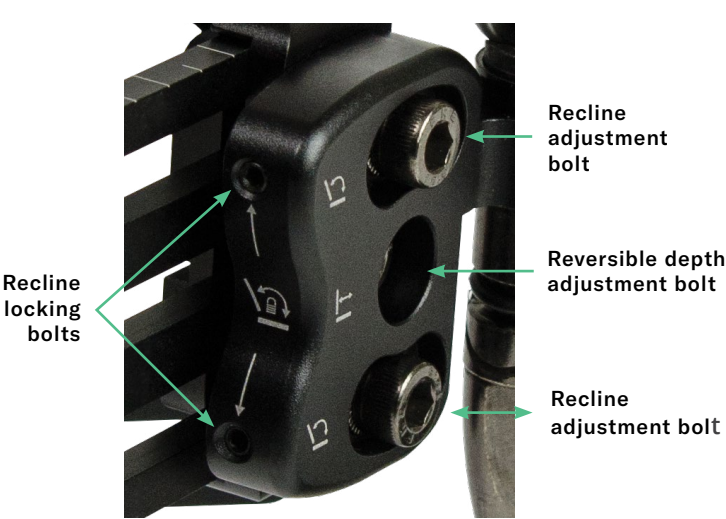

Reversible depth adjustment bolt allows for easy access even in the flipped cane clamp configurations.

Hardware options

- Compass 4 Quick Release 2 Point
- Compass 4 Fixed 2 Point
- Compass 4 Quick Release 4 Point
- Compass 4 Fixed 4 Point

Compass Power Mount Bracket


- Easy installation with 8" height adjustability
- Integrates the Comfort Unite, Embrace and Relief back supports when used with power seating systems:
 - Permobil Corpus® 3G (PMB-PERMOBIL)
 - Quantum TRU-Balance 3 or 4 (PMB-QUANTUM)
 - Quickie 3.7 (PMB-QUICKIE)
 - ROVI with Motion Concepts (PMB-ROVI)
 - AVID Rehab (PMB-AVID)



Compass Power Mount Bracket

Permobil Direct Mount Hardware

- Seamlessly compatible with Permobil Corpus 3G seating system
- Easy installation with less than 1/4" loss of seat depth
- Part number (BS-PDMXXW)

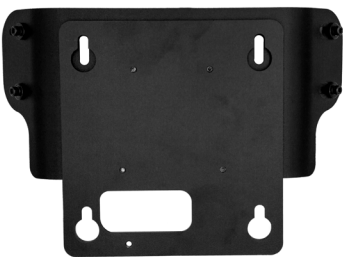

Permobil Direct Mount Hardware

Back support technology



The BOA® Fit System

- Micro-adjustable**
Each click of the BOA® dial offers one millimeter of precision fit adjustment.
- Secure**
Once the BOA® dial is locked into place, it stays in place.
- Customizable fit**
BOA® Technology delivers smooth, even adjustment with no pressure points for maximum comfort.
- Convenient**
The BOA® Fit System offers an adjustment solution free of messy straps, ratchets, and fasteners.

The BOA® dial




Pull to release tension



Push to engage the system
Turn right to tighten


Moldable support

As a standard offering on all Unite and Embrace back supports, the aluminum stays are located inside the zippered cover, behind the foam cushion surface and are accessible for adjustments while the individual is seated.

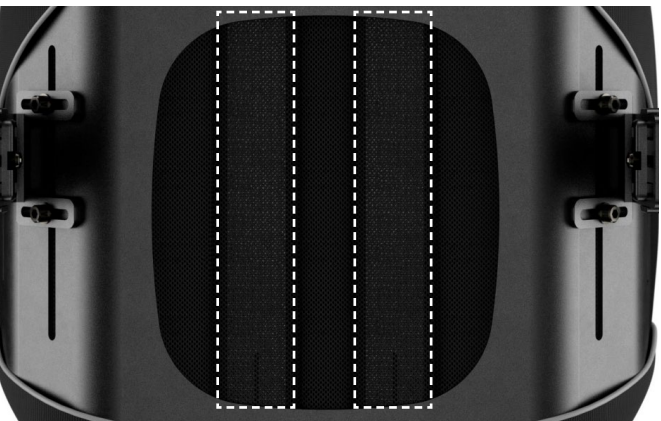

ROHO® air inserts

Available on Unite 1 Wing, Unite 2 and Unite 3 back supports, these adjustable air cell inserts are located along the spinal area and/or laterals for additional skin protection and positioning.

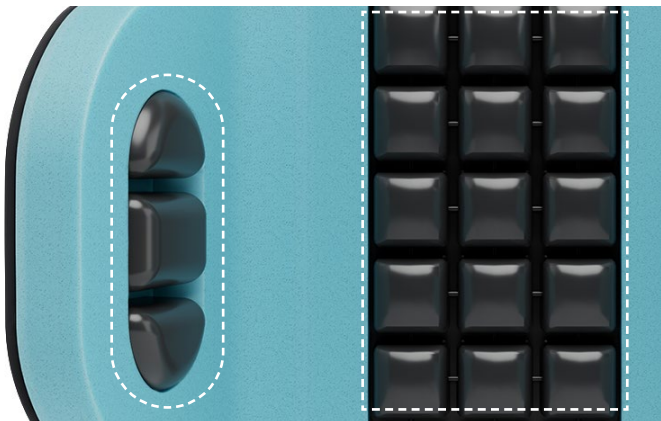

Custom back supports

Custom configurations available on the Comfort Unite, Relief and Embrace. Choose from foams, foam-in-place, ROHO air cells and additional size offerings.

E2617
Custom fabricated wheelchair back

Comfort Elements

E2611/E2612
General use back



Features

- Easily folds with wheelchair for transportation and storage
- Adjustable tension strapping system can accommodate for kyphosis and provide lumbar support
- 1” foam provides a comfortable back support
- Back securely installs directly over wheelchair canes
- Steel D-Rings provide superior strength and durability
- Cover options: Comfort-Tek or Stretch Air
- Warranty: 5 years back and hardware, 6 months for cover

Dimensions		Weight capacity	
Width	15” – 21”, 22” – 25”	All sizes	300 lbs
Length	16”		

Comfort Elements part number
ELB-BSXXXXW16L
<i>*Length (L) refers to the actual size dimension from bottom to top edge</i>
XXXX to be replaced XXXX = width in inches

Comfort Unite 1

E2613/E2614
Positioning back



Features

- 2.25" of shell contour.
- Adjustable moldable stays to conform and support the back.
- Keyhole (in the shell) and stud strap (on the cover) attachments for secure fit.
- Lightweight, durable aluminum shell providing pelvic and scapular relief.
- Hardware quick install: easy adjustment for depth, height and angle while user is in chair.
- Compatible with Comfort Compass hardware, Power Mount hardware and Direct Mount hardware.
- Cover options: Comfort-Tek or Diamond Spacer; GlideWear upgrade available on Diamond Spacer
- Available as a custom back through Comfort Designs.
- Warranty: 5 years back and hardware, 6 months for cover

Dimensions		Weight capacity	
Width	12” – 26”, (2" increments)	≤ 18” wide	250 lbs
Length	12” – 22”, (2" increments)	20” wide	300 lbs
Contours	2.25”	22” wide	350 lbs
		24” wide	400 lbs
		26" – 28" wide	450 lbs

Comfort Unite 1 part number
U1-BSXXWYYL
<i>*Length (L) refers to the actual size dimension from bottom to top edge</i>
XX & YY to be replaced XX = width & YY = length in inches
*See page 3 for hardware and cover options.

Standard sizes								
	12” W	14” W	16” W	18” W	20” W	22” W	24” W	26” W
22” L			16 x 22	18 x 22	20 x 22	22 x 22	24 x 22	
20” L		14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20	26 x 20
18” L	12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18	26 x 18
16” L	12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16		
14” L	12 x 14	14 x 14	16 x 14	18 x 14	20 x 14			
12” L	12 x 12	14 x 12	16 x 12	18 x 12				

Comfort Visco Back

E2611/E2612
General use back



Features

- Visco foam to redistribute pressure and create a supportive fit
- Three layer foam construction
- Lightweight and durable aluminum shell
- Compatible with head support and laterals
- Comfort Compass® hardware allows for easy height, depth and angle adjustment
- Cover options: Comfort-Tek or Stretch-Air
- Warranty: 5 years back and hardware, 6 months for cover

Dimensions		Weight capacity	
Width	16” – 24”, (2" increments)	≤ 18” wide	250 lbs
Length	14” – 22”, (2" increments)	20” wide	300 lbs
		22” wide	350 lbs
		24” wide	400 lbs

Comfort Visco part number
VB-BSXXWYYL
<i>*Length (L) refers to the actual size dimension from bottom to top edge</i>
XX & YY to be replaced XX = width & YY = length in inches
*See page 3 for hardware and cover options.

Standard sizes					
	16” W	18” W	20” W	22” W	24” W
22” L	16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20” L	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18” L	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16” L	16 x 16	18 x 16	20 x 16	22 x 16	
14” L	16 x 14	18 x 14	20 x 14		

Comfort Unite 2

E2615/E2616
Positioning back, Posterior



Dimensions		Weight capacity	
Width	12" – 24", (2" increments)	≤ 18" wide	250 lbs
Length	12" – 22", (2" increments)	20" wide	300 lbs
Contours	4"	22" wide	350 lbs
		24" wide	400 lbs

Standard sizes							
	12" W	14" W	16" W	18" W	20" W	22" W	24" W
22" L			16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L		14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14" L	12 x 14	14 x 14	16 x 14	18 x 14	20 x 14		
12" L	12 x 12	14 x 12	16 x 12	18 x 12			

Standard sizes with ROHO insert						
	14" W	16" W	18" W	20" W	22" W	24" W
22" L		16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L	14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14" L	14 x 14	16 x 14	18 x 14	20 x 14		

Features

- 4" of shell contour.
- Adjustable moldable stays to conform and support the back.
- Keyhole (in the shell) and stud strap (on the cover) attachments for secure fit.
- Lightweight, durable aluminum shell providing pelvic and scapular relief.
- Hardware quick install: easy adjustment for depth, height and angle while user is in chair.
- ROHO insert option along the spinal area.
- Compatible with Comfort Compass hardware, Power Mount hardware and Direct Mount hardware
- Cover options: Comfort-Tek or Diamond Spacer; GlideWear upgrade available on Diamond Spacer
- Available as a custom back through Comfort Designs.
- Warranty: 5 years back and hardware, 6 months for cover

Comfort Unite 2 part number	
U2-BSXXWYYL	
U2R-BSXXWYYL	With ROHO insert
<i>*Length (L) refers to the actual size dimension from bottom to top edge</i>	
XX & YY to be replaced XX = width & YY = length in inches	
<i>*See page 3 for hardware and cover options.</i>	

Comfort Unite 3

E2620/E2621
Positioning back support, Planar with lateral supports



Dimensions		Weight capacity	
Width	14" – 24", (2" increments)	≤ 18" wide	250 lbs
Length	14" – 22", (2" increments)	20" wide	300 lbs
Contours	6.5"	22" wide	350 lbs
		24" wide	400 lbs

Standard sizes						
	14" W	16" W	18" W	20" W	22" W	24" W
22" L		16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L	14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14" L	14 x 14	16 x 14	18 x 14	20 x 14		

Standard sizes with ROHO insert					
	16" W	18" W	20" W	22" W	24" W
22" L	16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	16 x 16	18 x 16	20 x 16	22 x 16	

Features

- 6.5" of shell contour.
- Adjustable moldable stays to conform and support the back.
- Keyhole (in the shell) and stud strap (on the cover) attachments for secure fit.
- Lightweight, durable aluminum shell providing pelvic and scapular relief.
- ROHO insert option along the spine and in laterals.
- Hardware quick install: easy adjustment for depth, height and angle while user is in chair.
- Compatible with Comfort Compass hardware, Power Mount hardware and Direct Mount hardware.
- Cover options: Comfort-Tek or Diamond Spacer; GlideWear upgrade available on Diamond Spacer
- Available as a custom back through Comfort Designs.
- Warranty: 5 years back and hardware, 6 months for cover

Comfort Unite 3 part number	
U3-BSXXWYYL	
U3R-BSXXWYYL	With ROHO insert
<i>*Length (L) refers to the actual size dimension from bottom to top edge</i>	
XX & YY to be replaced XX = width & YY = length in inches	
<i>*See page 3 for hardware and cover options.</i>	

Comfort Unite 1 Wing

E2620/E2621
Positioning back, Posterior



Features

- 2.25" of shell contour, available with 5" or 7" deep wing.
- Adjustable moldable stays to conform and support the back.
- Keyhole (in the shell) and stud strap (on the cover) attachments for secure fit.
- Lightweight, durable aluminum shell providing pelvic and scapular relief.
- ROHO insert option along the spinal area.
- Hardware quick install: easy adjustment for depth, height and angle while user is in chair.
- Compatible with Comfort Compass hardware, Power Mount hardware and Direct Mount hardware.
- Cover options: Comfort-Tek or Diamond Spacer; GlideWear upgrade available on Diamond Spacer
- Available as a custom back through Comfort Designs.
- Warranty: 5 years back and hardware, 6 months for cover

Dimensions		Weight capacity	
Width	12" – 24", (2" increments)	≤ 18" wide	250 lbs
Length	12" – 22", (2" increments)	20" wide	300 lbs
Contours	2.25" with 5" or 7" wing	22" wide	350 lbs
		24" wide	400 lbs

Unite 1 Wing 5" standard sizes

	12" W	14" W	16" W	18" W	20" W	22" W	24" W
22" L			16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L		14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14" L	12 x 14	14 x 14	16 x 14	18 x 14	20 x 14		
12" L	12 x 12	14 x 12	16 x 12	18 x 12			

Unite 1 Wing 7" standard sizes

	12" W	14" W	16" W	18" W	20" W	22" W	24" W
22" L			16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L		14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14" L			16 x 14	18 x 14	20 x 14		
12" L			16 x 12	18 x 12			

Comfort Unite 1 Wing part number

U1W5-BSXXWYYL (5" WING)

U1W7-BSXXWYYL (7" WING)

U1W5R-BSXXWYYL (5" WING) With ROHO insert

U1W7R-BSXXWYYL (7" WING) With ROHO insert

**Length (L) refers to the actual size dimension from bottom to top edge*

XX & YY to be replaced | XX = width & YY = length in inches

**See page 3 for hardware and cover options.*

Unite 1 Wing 5" standard sizes with ROHO insert

	14" W	16" W	18" W	20" W	22" W	24" W
22" L		16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L	14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14" L	14 x 14	16 x 14	18 x 14	20 x 14		

Unite 1 Wing 7" standard sizes with ROHO insert

	14" W	16" W	18" W	20" W	22" W	24" W
22" L		16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20" L	14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18" L	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16" L	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14" L		16 x 14	18 x 14	20 x 14		

Comfort Relief

E2615
Positioning back, Posterior lateral



Features

- The open back design gives the ability to off-load the pressure typically seen on spinal deformities through the back shell opening.
- With the BOA® Fit System, equal and stable pressure distribution is accomplished in an “off the shelf” back design.
- Updated cover design and 3-layer foam construction.
- BodiLink® and hardware are optional and available.
- Cover options: Comfort-Tek and Diamond Spacer
- Available as a custom back through Comfort Designs.
- Warranty: 5 years back and hardware, 6 months for cover

Comfort Relief part number

R-BSXXWYYL

**Length (L) refers to the actual size dimension from bottom to top edge*

XX & YY to be replaced | XX = width & YY = length in inches

**See page 3 for hardware and cover options.*

Comfort Relief head support bracket

Back support width	Part number
14W	BL-HS-AR14W
16W	BL-HS-AR16W
18W	BL-HS-AR18W
20W	BL-HS-AR20W


Standard sizes

	14" W	16" W	18" W	20" W
22" L	14 x 22	16 x 22	18 x 22	20 x 22
20" L	14 x 20	16 x 20	18 x 20	20 x 20
18" L	14 x 18	16 x 18	18 x 18	20 x 18
16" L	14 x 16	16 x 16	18 x 16	20 x 16
14" L	14 x 14	16 x 14	18 x 14	20 x 14
12" L	14 x 12	16 x 12	18 x 12	20 x 12
10" L	14 x 10	16 x 10	18 x 10	20 x 10

THE BOA® dial

See page 6 for details.

Pressure redistribution properties

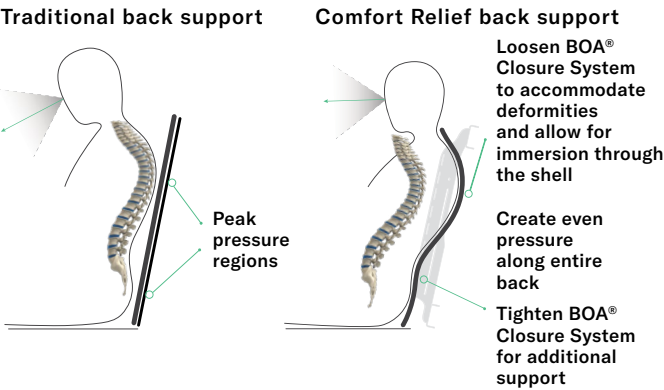

Comfort Embrace

E2613 / E2614
Positioning back, Posterior



Features

- Integrated wings attach to a planar aluminum back support shell and are adjustable for angle, height, and up to a total of 3” of width for standard wings and 1.5” for pediatric wings.
- Moldable stays included for positioning, support and pressure redistribution for a variety of spinal curvatures and changes in conditions.
- Updated cover design and 2 layer foam construction.
- Comes with integrated chest harness strap guides.
- Durable aluminum shell comes with pre-drilled 2” x 1” head support hole pattern for BodiLink Head Support.
- Compatible with lateral trunk support hardware. A total of 11” of pad width adjustment (5.5” per side) & over 5” of height adjustment.
- Cover options: Comfort-Tek and Diamond Spacer
- Available as a custom back through Comfort Designs.
- Warranty: 5 years back and hardware, 6 months for cover

Dimensions		Weight capacity	
Width	10” – 24”, (2" increments)	10” - 18” wide	250 lbs
Length	10” – 22”, (2" increments)	20” wide	300 lbs
		22” wide	350 lbs
		24” wide	400 lbs

Standard sizes								
	10” W	12” W	14” W	16” W	18” W	20” W	22” W	24” W
22” L				16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20” L			14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18” L		12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16” L		12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14” L	10 x 14	12 x 14	14 x 14	16 x 14	18 x 14	20 x 14		
12” L	10 x 12	12 x 12	14 x 12	16 x 12	18 x 12			
10” L	10 x 10	12 x 10	14 x 10	16 x 10				

*Length (L) refers to the actual size dimension of the support from bottom to top edge.
*24” Wide back supports require 4 Point hardware or Power Mount hardware.
*Power Mount hardware not available on 10”W or 10”L back supports.
*4 Point hardware not available on 10” L if using head support.

ROHO Agility Active Aluminum

E2620/E2621
Positioning back, planar with lateral supports



Dimensions		Weight capacity	
Width	14" - 18" (1" increments)	All sizes	265 lbs
Length	8", 10"		
Contours	Min: 3" ; Mid: 4"		

Standard sizes					
	14” W	15” W	16” W	17” W	18” W
10” L	10x14	10x15	10x16	10x17	10x18
8” L	8x14	8x15	8x16	8x17	8x18

Features

- Adjustable hardware includes tough, durable Zytel® discs that ratchet and lock into the aluminum shell.
- The adjustable discs rotate independently providing 15° of forward and rearward angle adjustment in 5° increments.
- The shell can be flipped to accommodate individual body types.
- Warranty: 5 years back shell & hardware, 6 months cover warranty

ROHO Agility Active Aluminum part number

AG2- WWLLF -2A	Minimum contour
AG3- WWLLF -2A	Mid contour

*Length (L) refers to the actual size dimension from bottom to top edge

WW & **LL** to be replaced | **WW** = width & **LL** = length in inches

ROHO Agility Active Carbon

E2620/E2621
Positioning back, planar with lateral supports



Dimensions		Weight capacity	
Width	14" - 18" (1" increments)	All sizes	265 lbs
Length	10", 13", 16"		
Contours	Min: 3" ; Mid: 4"		

Standard sizes					
	14” W	15” W	16” W	17” W	18” W
16” L	16x14	16x15	16x16	16x17	16x18
13” L	13x14	13x15	13x16	13x17	13x18
10” L	10x14	10x15	10x16	10x17	10x18

Features

- Flat-weave carbon fiber composition lays closer together to increase its strength while remaining lightweight.
- Adjustable hardware includes tough, durable Zytel discs that ratchet and lock into the aluminum shell.
- The adjustable discs rotate independently providing 15° of forward and rearward angle adjustment in 5° increments.
- Includes adjustable air insert along spinal area.
- The shell can be flipped to accommodate individual body types.
- Warranty: 5 years back shell & hardware, 6 months cover warranty

ROHO Agility Active Carbon part number

AG2- WWLLA -2CB	Minimum contour
AG3- WWLLA -2CB	Mid contour

*Length (L) refers to the actual size dimension from bottom to top edge

WW & **LL** to be replaced | **WW** = width & **LL** = length in inches

Back supports with included trunk supports

- Lateral trunk supports (LTS) are constructed with moldable stainless steel and protective foam.
- Maximum distance between two lateral supports is equal to the width of the wheelchair. Minimum distance between two lateral supports is 8” less than the width of the wheelchair.
- Lateral trunk supports have up to 8” of width adjustment (4” per pad) and 1.5” depth adjustment each.
- Lateral trunk supports are available in four different sizes: extra small, small, medium, and large.
- Fixed lateral hardware comes standard.
- Extra lateral hardware links available to provide optimal contact with the user.

Lateral trunk support options

- 4” D x 3” L (LAT-4D3L-COMFORT-TEK or LAT-4D3L-STRETCH-AIR)
- 5” D x 4” L (LAT-5D4L-COMFORT-TEK or LAT-5D4L-STRETCH-AIR)
- 6” D x 4” L (LAT-6D4L-COMFORT-TEK or LAT-6D4L-STRETCH-AIR)
- 7” D x 5” L (LAT-7D5L-COMFORT-TEK or LAT-7D5L-STRETCH-AIR)

Lateral trunk support option can be selected by adding the (part number) to the Comfort Unite 1 LTS, Comfort Relief LTS or Comfort Embrace LTS part number.

Comfort Unite 1 LTS



Weight capacity	
10” - 18” wide	250 lbs
20” wide	300 lbs
22” wide	350 lbs
24” wide	400 lbs

E2620/E2621
Positioning back,
Planar with lateral supports

Comfort Unite LTS part number	
U1LTS-BSXXWYYL	
*Length (L) refers to the actual size dimension from bottom to top edge	
XX & YY to be replaced XX = width & YY = length in inches	
*See page 3 for hardware and cover options.	

Standard sizes							
	12” W	14” W	16” W	18” W	20” W	22” W	24” W
22” L			16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20” L		14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18” L		14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16” L		14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14” L	12 x 14	14 x 14	16 x 14	18 x 14	20 x 14		
12” L	12 x 12	14 x 12	16 x 12	18 x 12			

Comfort Relief LTS



Weight capacity	
14” - 18” wide	250 lbs
20” wide	300 lbs

Comfort Embrace LTS



Weight capacity	
10” - 18” wide	250 lbs
20” wide	300 lbs
22” wide	350 lbs
24” wide	400 lbs

E2620
Positioning back,
Planar with lateral supports

Comfort Relief LTS part number	
RLTS-BSXXWYYL	
*Length (L) refers to the actual size dimension from bottom to top edge	
XX & YY to be replaced XX = width & YY = length in inches	
*See page 3 for hardware and cover options.	

Comfort Relief head support bracket	
Back support width	Part number
14W	BL-HS-AR14W
16W	BL-HS-AR16W
18W	BL-HS-AR18W
20W	BL-HS-AR20W

Standard sizes				
	14” W	16” W	18” W	20” W
22” L	14 x 22	16 x 22	18 x 22	20 x 22
20” L	14 x 20	16 x 20	18 x 20	20 x 20
18” L	14 x 18	16 x 18	18 x 18	20 x 18
16” L	14 x 16	16 x 16	18 x 16	20 x 16
14” L	14 x 14	16 x 14	18 x 14	20 x 14
12” L	14 x 12	16 x 12	18 x 12	20 x 12

E2620 / E2621
Positioning back,
Planar with lateral supports

Comfort Embrace LTS part number	
ELTS-BSXXWYYL	
*Length (L) refers to the actual size dimension from bottom to top edge	
XX & YY to be replaced XX = width & YY = length in inches	
*See page 3 for hardware and cover options.	

Standard sizes								
	10” W	12” W	14” W	16” W	18” W	20” W	22” W	24” W
22” L				16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
20” L			14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
18” L		12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
16” L		12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	
14” L	10 x 14	12 x 14	14 x 14	16 x 14	18 x 14	20 x 14		
12” L	10 x 12	12 x 12	14 x 12	16 x 12	18 x 12			
10” L	10 x 10	12 x 10	14 x 10	16 x 10				



Cushions

Permobil cushions are engineered and designed to provide superior pressure relief and positioning to accommodate a wide range of seating requirements. The utilization of several patented technologies allows users to experience an unmatched level of comfort and independence designed to enhance their everyday lives.

GlideWear features

Dual-layer Shear Reduction Zone™

Greatly reducing the effects of shear forces, a known contributor to the breakdown of tissue.

Micro-climate control

GlideWear is lightweight, breathable and reduces moisture buildup for micro-climate control.

Stability

Located under the bony prominences, it allows for micro-movements without the deformation of tissue in the most at-risk areas.

Works in combination with the higher friction areas of the seating surface surrounding it, keeping the individual stable.




Seating dimensions

- A – Seat thickness

B – Inferior thigh support thickness
- C – Medial thigh support height


D – Lateral thigh support height

E – Pelvic contour height

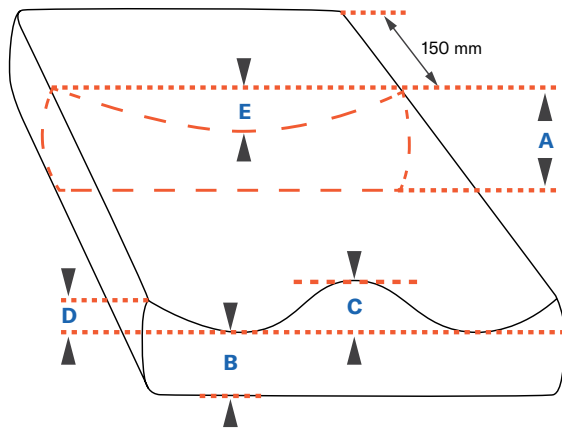

QuadraGel®

The pressure-relieving QuadraGel pack reduces heat build up. Four compartments evenly distribute gel under the seating surface.

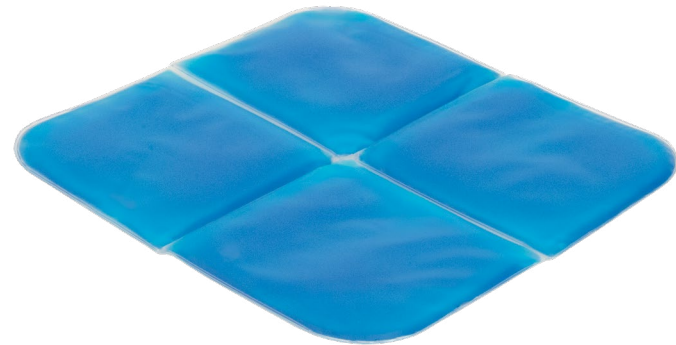

Adjustable Technology Insert (ATI)

Constructed of a strong aluminum that is designed to be molded to add additional support where needed. Allows for on-site adjustments to effectively meet the needs of the individual.

The insert can be remolded to meet the changes in the individual's physical and functional needs.



**See instruction manual for proper utilization*

ROHO Mosaic



Dimensions		Weight capacity	
Cell height	2.5 – 3"	All sizes	315 lbs

Standard sizes	16" W	17" W	18" W	19" W	20" W
16" D	16x16		18x16		
17" D		17x17			
18" D	16x18		18x18		20x18

Elements with gel



Weight capacity	
16" wide	250 lbs
18" wide	300 lbs
20" wide	350 lbs
22" wide & wider	400 lbs

Standard sizes	16" W	18" W	20" W	22" W	24" W
16" D	16 x 16	18 x 16	20 x 16	22 x 16	24 x 16
18" D	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
20" D	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20

E2601
General use

Features

- Interconnected large air cells.
- Available in six sizes.
- Cover options: Standard or Heavy Duty.
- Warranty: 6 months cover / 1 year cushion

Mosaic part numbers

Cover style	Standard cover	Heavy duty cover
	MOS WWDD CA	MOS WWDD DH

WWDD to be replaced | **WW** = width & **DD** = depth in inches

QuadraGel®

The QuadraGel prevents gel migration and redistributes pressure for comfort.



E2601/E2602
General use

Features

- Non-skid flat bottom fabric with available hook and loop attachments ([HOOK-KIT-4](#)).
- Single layer foam construction contours to the user's body reducing high pressure areas.
- Warranty: 6 months cover / 5 years cushion

Elements with gel part numbers

Seat thickness	Stretch-Air	Comfot-Tek
2" Seat thickness	46GS- WWDD -B	46G- WWDD -B
3" Seat thickness	463GS- WWDD -B	463G- WWDD -B

WWDD to be replaced | **WW** = width & **DD** = depth in inches

Comfort Support Pro 7 Series

E2602
General use



Features

- Non-contoured molded foam surface allowing for easier transfers.
- Coccyx relief to assist in offloading bony prominences.
- Kwik strap loops around wheelchair canes to securely attach cushion during transfers.
- Warranty: 6 months cover/ 5 years cushion

Dimensions		Weight capacity	
Seat thickness	3.5"	All sizes	750 lbs

Standard sizes					
	22" W	24" W	26" W	28" W	30" W
18" D	22 x 18	24 x 18	26 x 18	28 x 18	30 x 18
20" D	22 x 18	24 x 20	26 x 20	28 x 20	30 x 20

Support Pro 7 Series part number	
53- WWDD -B	Without QuadraGel
53G- WWDD -B	With QuadraGel
WW & DD to be replaced WW = width & DD = depth in inches	

Comfort Ascent

E2603 / E2604
Skin protection



Features

- Tapered leg adductors increase in width near the back of the cushion to give support to the greater trochanters.
- Layered high resiliency foams, with an EVA foam bottom layer, create a soft and durable seating surface.
- Warranty: 6 months cover/ 5 years cushion

Ascent part numbers	
Cover style	Comf-Tek
Stretch-Air	AC-GF- WWDD
WWDD to be replaced WW = width & DD = depth in inches	

Standard sizes					
	16" W	18" W	20" W	22" W	24" W
16" D	16 x 16	18 x 16	20 x 16		
18" D	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
20" D	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20

*Odd sizes available in 1-inch increments.

Curve®

E2601/E2602
General use



Dimensions		Weight capacity	
Seat thickness	3"	12" wide & under	100 lbs
Inferior thigh support thickness	2.25"	14" wide	150 lbs
Medial thigh support thickness	1"	16" wide	250 lbs
Lateral thigh support height	1"	18" wide	300 lbs
Pelvic contour height	1"	20" wide	350 lbs
		22" wide & wider	400 lbs

Standard sizes								
	8" W	10" W	12" W	14" W	16" W	18" W	20" W	24" W
8" D	8 x 8							
10" D	8 x 10	10 x 10	12 x 10	14 x 10				
12" D	8 x 12	10 x 12	12 x 12	14 x 12				
14" D		10 x 14	12 x 14	14 x 14	16 x 14			
16" D			12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	24 x 16
18" D				14 x 18	16 x 18	18 x 18	20 x 18	24 x 18
20" D					16 x 20	18 x 20	20 x 20	24 x 20

*Odd sizes available in 1-inch increments.

Curve part numbers		
Cover style	Stretch-Air	Comf-Tek
Standard	CU-SV- WWDD	CU-FV- WWDD
WWDD to be replaced WW = width & DD = depth in inches		

ROHO AirLite®

E2605
Positioning



Features

- Cover is fluid-resistant, machine washable and can be disinfected.
- Lower anterior height for ease with transfers.
- Non-adjustable ROHO Air Floatation air insert built into the cushion.
- Warranty: 6 months cover / 2 years cushion

Standard sizes							
	13" W	14" W	15" W	16" W	17" W	18" W	20" W
13" D	13 x 13						
14" D		14 x 14					
15" D	13 x 15		15 x 15				
16" D		14 x 16		16 x 16		18 x 16	
17" D			15 x 17		17 x 17		20 x 17
18" D				16 x 18	17 x 18	18 x 18	20 x 18
19" D							19 x 19
20" D					17 x 20	18 x 20	

AirLite part number	
AL WWDD HD	
WW & DD to be replaced WW = width & DD = depth in inches	

Dimensions		Weight capacity	
Cell heights	2" - 3.75"	All sizes	300 lbs
Contoured foam	Yes		

Comfort Saddle 7 Series



Dimensions		Weight capacity	
Seat thickness	4.5"	All sizes	750 lbs
Inferior thigh support thickness	3.5"		
Abductor	5"		

Standard sizes	22" W	24" W	26" W	28" W	30" W
18" D	22 x 18	24 x 18	26 x 18	28 x 18	30 x 18
20" D	22 x 18	24 x 20	26 x 20	28 x 20	30 x 20

E2608
Skin protection and positioning

Features

- Lateral and medial contours promote postural alignment and stability.
- Optional QuadraGel pack to redistribute pressure for skin protection.
- Coccyx relief to assist in offloading bony prominences.
- Kwik strap loops around wheelchair canes to securely attach cushion during transfers.
- Warranty: 6 months cover/ 5 years cushion

Saddle 7 Series part number	
53S- WWDD -B	Without QuadraGel
53SG- WWDD -B	With QuadraGel
WW & DD to be replaced WW = width & DD = length in inches	

ROHO Hybrid Elite SR® Single Compartment



E2622 / E2623
Adjustable skin protection

Features

- Pelvic loading area is adjustable to shape, weight, and changes in individual's position over time.
- The air-filled cushion insert can be easily cleaned and disinfected.
- Comes Sensor Ready for Smart Check® compatibility
- Warranty: 6 months cover / 2 years cushion

Weight capacity	Dimensions
All sizes	500 lbs
	Cell height
	4.25"

Standard sizes	14" W	15" W	16" W	17" W	18" W	19" W	20" W	21" W	22" W	23" W	24" W
14" D	14 x 14										
15" D	14 x 15	15 x 15	16 x 15								
16" D	14 x 16	15 x 16	16 x 16	17 x 16	18 x 16						
17" D		15 x 17	16 x 17	17 x 17	18 x 17		20 x 17				
18" D			16 x 18	17 x 18	18 x 18		20 x 18		22 x 18		24 x 18
19" D											
20" D			16 x 20	17 x 20	18 x 20		20 x 20		22 x 20		24 x 20
21" D											
22" D									22 x 22		

Smart Check part number
SMARTCHECK-01

ROHO Hybrid Elite Dual Compartment



Weight capacity	Dimensions
All sizes	500 lbs
	Cell height
	4.25"

Standard sizes	14" W	15" W	16" W	17" W	18" W	19" W	20" W	21" W	22" W	23" W	24" W
14" D	14 x 14										
15" D		15 x 15	16 x 15								
16" D	14 x 16	15 x 16	16 x 16	17 x 16	18 x 16						
17" D		15 x 17	16 x 17	17 x 17	18 x 17		20 x 17				
18" D			16 x 18	17 x 18	18 x 18		20 x 18		22 x 18		24 x 18
19" D											
20" D			16 x 20	17 x 20	18 x 20		20 x 20		22 x 20		24 x 20
21" D											
22" D									22 x 22		

Features

- Pelvic loading area, with dual compartments, is adjustable to shape, weight, and changes in individual's position over time.
- Contoured foam base and air insert are fluid resistant.
- Air-filled cushion insert can be easily cleaned and disinfected.
- Warranty: 6 months cover / 2 years cushion

Hybrid Elite Dual Compartment part numbers		
Standard cover	Heavy duty cover >20"	Heavy duty cover <22"
2RHE WWDD C	2RHE WWDD DHD	2RHE WWDD HDLG
WWDD to be replaced WW = width & DD = depth in inches		

ROHO Single Compartment

E2622 / E2623
Adjustable skin protection



Single Compartment part numbers

Cushion style	Standard cover	Heavy duty cover
High profile	1R WWDDC	1R WWDDHD
Mid profile	1R WWDDMPC	-
Low profile	1R WWDDLPC	-

	Sensor ready standard cover	Sensor ready heavy duty cover
High profile	1R WWDDH -CA-SR	1R WWDDH -CH-SR
Mid profile	1R WWDDM -CA-SR	R WWDDM -CH-SR

WWDD to be replaced | **WW** = width & **DD** = depth in cells

Features

- Adjustable to shape, weight, and changes in an individual’s position over time.
- Hydrostatic distribution of load across the entire contact area.
- Provides immersion and envelopment for pressure redistribution.
- SR: Cushion is Sensor Ready™ for Smart Check, which can be easily attached and detached by means of quick disconnect.
- No weight limit if appropriately sized to the individual.
- Warranty: 6 months cover / 2 years cushion

Smart Check part number

SMARTCHECK-01

Dimensions

High profile cell height	4.25"
Mid profile cell height	3.25"
Low profile cell height	2.25"

High profile standard cover sizes

	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W	14 W	15 W
6 D	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6	14 x 6	15 x 6
7 D	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7	14 x 7	15 x 7
8 D	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8	14 x 8	15 x 8
9 D	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9	14 x 9	15 x 9
10 D	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10	14 x 10	15 x 10
11 D	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11	14 x 11	15 x 11
12 D	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12	14 x 12	15 x 12
13 D	6 x 13	7 x 13	8 x 13	9 x 13	10 x 13	11 x 13	12 x 13	13 x 13	14 x 13	15 x 13
14 D	6 x 14	7 x 14	8 x 14	9 x 14	10 x 14	11 x 14	12 x 14	13 x 14	14 x 14	15 x 14

High profile heavy duty cover sizes

	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W	14 W	15 W
6 D	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6	14 x 6	15 x 6
7 D	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7	14 x 7	15 x 7
8 D	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8	14 x 8	15 x 8
9 D	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9	14 x 9	15 x 9
10 D	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10	14 x 10	15 x 10
11 D	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11	14 x 11	15 x 11
12 D	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12	14 x 12	15 x 12
13 D	6 x 13	7 x 13	8 x 13	9 x 13	10 x 13	11 x 13	12 x 13	13 x 13	14 x 13	15 x 13
14 D	6 x 14	7 x 14	8 x 14	9 x 14	10 x 14	11 x 14	12 x 14	13 x 14	14 x 14	15 x 14

Sizing key

	Standard sizes
	Made to order

High profile sensor ready sizes

	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W	14 W	15 W
6 D	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6	14 x 6	15 x 6
7 D	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7	14 x 7	15 x 7
8 D	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8	14 x 8	15 x 8
9 D	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9	14 x 9	15 x 9
10 D	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10	14 x 10	15 x 10
11 D	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11	14 x 11	15 x 11
12 D	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12	14 x 12	15 x 12

Mid profile sizes

	7 W	8 W	9 W	10 W	11 W	12 W	13 W
7 D	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7
8 D	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8
9 D	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9
10 D	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10
11 D	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11
12 D	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12

Mid profile sensor ready sizes

	7 W	8 W	9 W	10 W	11 W	12 W	13 W
7 D	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7
8 D	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8
9 D	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9
10 D	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10
11 D	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11
12 D	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12

Low profile sizes

	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W
6 D	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6
7 D	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7
8 D	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8
9 D	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9
10 D	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10
11 D	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11
12 D	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12

Note: These cushion sizes are not in inches, but in cells. Please view the cell conversion table on page 34, to see the sizing of cells to inches.

ROHO Hybrid Select



Weight capacity		Dimensions	
All sizes	500 lbs	Cell height	2.5" - 5.13" contoured

Standard sizes							
	14" W	15" W	16" W	17" W	18" W	19" W	20" W
14" D	14 x 14						
15" D		15 x 15	16 x 15				
16" D	14 x 16	15 x 16	16 x 16	17 x 16	18 x 16		20 x 16
17" D		15 x 17	16 x 17	17 x 17	18 x 17		20 x 17
18" D			16 x 18	17 x 18	18 x 18		20 x 18
19" D							
20" D			16 x 20	17 x 20	18 x 20		20 x 20

ROHO Dual Compartment



Dimensions	
High profile cell height	4.25"
Low profile cell height	2.25"

ROHO sizing key	
	Standard sizes
	Made to order

Note: These cushion sizes are not in inches, but in cells. Please view the cell conversion table on page 34, to see the sizing of cells to inches.

High profile sizes				
	8 W	9 W	10 W	11 W
8 D	8 x 8	9 x 8	10 x 8	11 x 8
9 D	8 x 9	9 x 9	10 x 9	11 x 9
10 D	8 x 10	9 x 10	10 x 10	11 x 10
11 D	8 x 11	9 x 11	10 x 11	11 x 11
12 D	8 x 12	9 x 12	10 x 12	11 x 12

E2624
Adjustable skin protection and positioning

Features

- Provides stability and positioning for the pelvis and thighs with a unique ROHO air cell design and contoured foam base
- Warranty: 6 months cover / 3 years cushion

Hybrid Select part numbers

Cushion style	Standard cover
High profile	HS WWDDC
WWDD to be replaced WW = width & DD = depth in inches	

ROHO Quadtro Select®



Dimensions	
High profile cell height	4.25"
Mid profile cell height	3.25"
Low profile cell height	2.25"

High profile standard cover sizes									
	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W	14 W
6 D	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6	14 x 6
7 D	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7	14 x 7
8 D	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8	14 x 8
9 D	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9	14 x 9
10 D	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10	14 x 10
11 D	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11	14 x 11
12 D	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12	14 x 12
13 D	6 x 13	7 x 13	8 x 13	9 x 13	10 x 13	11 x 13	12 x 13	13 x 13	14 x 13
14 D	6 x 14	7 x 14	8 x 14	9 x 14	10 x 14	11 x 14	12 x 14	13 x 14	14 x 14

High profile heavy duty cover sizes									
	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W	14 W
6 D	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6	14 x 6
7 D	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7	14 x 7
8 D	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8	14 x 8
9 D	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9	14 x 9
10 D	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10	14 x 10
11 D	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11	14 x 11
12 D	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12	14 x 12
13 D	6 x 13	7 x 13	8 x 13	9 x 13	10 x 13	11 x 13	12 x 13	13 x 13	14 x 13
14 D	6 x 14	7 x 14	8 x 14	9 x 14	10 x 14	11 x 14	12 x 14	13 x 14	14 x 14

Low profile standard cover sizes									
	5 W	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W
6 D	5 x 6	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6
7 D	5 x 7	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7
8 D	5 x 8	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8
9 D	5 x 9	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9
10 D	5 x 10	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10
11 D	5 x 11	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11
12 D	5 x 12	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12

Note: These cushion sizes are not in inches, but in cells. Please view the cell conversion table on page 34, to see the sizing of cells to inches.

E2624 / E2625
Adjustable skin protection and positioning

Features

- ISOFOLO Memory Control separates the cushion into quadrants which holds the pelvis left to right and front to back
- Adjustable to shape, weight, and changes in individual's position over time.
- Hydrostatic distribution of load across the entire contact area.
- No weight limit if appropriately sized to the individual.
- Warranty: 6 months cover / 3 years cushion

Quadtro Select part numbers

Cushion style	Standard cover	Heavy duty cover
High profile	QS WWDDC	QS WWDDHD
Mid profile	QS WWDDMPC	-
Low profile	QS WWDDLPC	-
WWDD to be replaced WW = width & DD = depth in cells		

Mid profile standard cover sizes

	7 W	8 W	9 W	10 W	11 W
7 D	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7
8 D	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8
9 D	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9
10 D	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10
11 D	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11

Mid profile heavy duty cover sizes

	7" W	8" W	9" W	10" W	11" W
7" D	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7
8" D	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8
9" D	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9
10" D	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10
11" D	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11

Low profile heavy duty cover sizes									
	5 W	6 W	7 W	8 W	9 W	10 W	11 W	12 W	13 W
6 D	5 x 6	6 x 6	7 x 6	8 x 6	9 x 6	10 x 6	11 x 6	12 x 6	13 x 6
7 D	5 x 7	6 x 7	7 x 7	8 x 7	9 x 7	10 x 7	11 x 7	12 x 7	13 x 7
8 D	5 x 8	6 x 8	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8
9 D	5 x 9	6 x 9	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9
10 D	5 x 10	6 x 10	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10
11 D	5 x 11	6 x 11	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11
12 D	5 x 12	6 x 12	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12

ROHO Contour Select®




Dimensions

Cell height	2.25" - 4.25" contoured
-------------	-------------------------

Contour Select sizes

	7 W	8 W	9 W	10 W	11 W	12 W	13 W
8 D	7 x 8	8 x 8	9 x 8	10 x 8	11 x 8	12 x 8	13 x 8
9 D	7 x 9	8 x 9	9 x 9	10 x 9	11 x 9	12 x 9	13 x 9
10 D	7 x 10	8 x 10	9 x 10	10 x 10	11 x 10	12 x 10	13 x 10
11 D	7 x 11	8 x 11	9 x 11	10 x 11	11 x 11	12 x 11	13 x 11
12 D	7 x 12	8 x 12	9 x 12	10 x 12	11 x 12	12 x 12	13 x 12

ROHO Enhancer

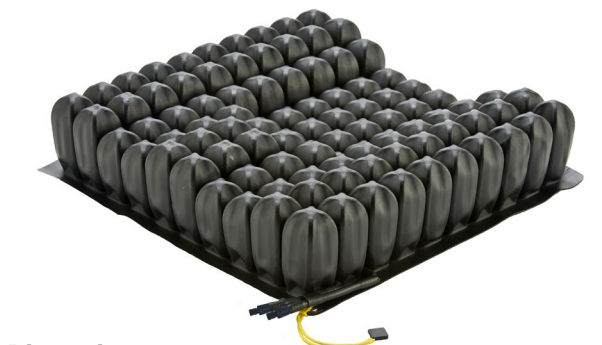

Dimensions

Cell height	2.25" - 4.25" contoured
-------------	-------------------------

Enhancer sizes

	8 W	9 W	10 W	11 W
8 D	8 x 8	9 x 8	10 x 8	11 x 8
9 D	8 x 9	9 x 9	10 x 9	11 x 9
10 D	8 x 10	9 x 10	10 x 10	11 x 10
11 D	8 x 11	9 x 11	10 x 11	11 x 11
12 D	8 x 12	9 x 12	10 x 12	11 x 12

ROHO sizing key

	Standard sizes
	Made to order

Note: These cushion sizes are not in inches, but in cells. Please view the cell conversion table on page 34, to see the sizing of cells to inches.

E2624 / E2625
Adjustable skin protection and positioning

Features

- The pre-contoured cells of the cushion assist in guiding the pelvis into the center of the ischial well and the thighs into the leg channels.
- Adjustable to shape, weight, and changes in an individual’s position over time.
- ISOFLO Memory Control separates the cushion into quadrants which holds the pelvis left to right and front to back.
- No weight limit if appropriately sized to the individual.
- Warranty: 6 months cover / 3 years cushion

Contour Select part numbers

Standard cover
CS WWDDC
WWDD to be replaced WW = width & DD = depth in cells

Features

- Adjustable to shape, weight, and changes in an individual’s position over time.
- Multiple cell heights & dual compartment system.
- Offers the ability to hold the pelvis and the thighs while offering adjustable lateral stability.
- Allows you to decrease the pressure under the ITs and enhance the lateral support along the greater trochanters.
- No weight limit if appropriately sized to the individual.
- Warranty: 6 months cover / 2 years cushion

Enhancer part numbers

Standard cover
EN HWDDC
WWDD to be replaced WW = width & DD = depth in cells

Comfort Embrace & Comfort Embrace ATI

E2607 / E2608
Skin protection & positioning
E2624 / E2625
Adjustable skin protection & positioning



ATI insert

Features

- Deep lateral and medial leg contours promote postural alignment.
- Tapered leg adductors increase in width near back of cushion to give support to the greater trochanters.
- Two foam layers provide comfort and skin protection: a contoured molded base foam and a layer of memory foam.
- Non-skid bottom fabric with available hook and loop attachments ([HOOK-KTT-4](#)).
- Warranty: 6 months cover / 5 years cushion

Dimensions	Zero Elevation	Anti-Thrust
16" x 14" and smaller		
Seat thickness	3"	2.75"
Inferior thigh support thickness	2.5"	3"
Medial thigh support height	.75"	.75"
Lateral thigh support height	1"	1"
Pelvic contour height	1"	1"

Dimensions	Zero Elevation	Anti-Thrust
16" x 16" and larger		
Seat thickness	4.25"	3.75"
Inferior thigh support thickness	3"	3.75"
Medial thigh support height	1.25"	1.25"
Lateral thigh Support height	1.5"	1.5"
Pelvic contour height	1.5"	1.25"

Available styles

Zero Elevation

Anti-Thrust



Standard sizes

	10" W	12" W	14" W	16" W	18" W	20" W	22" W	24" W
10" D	10 x 10	12 x 10	14 x 10					
12" D	10 x 12	12 x 12	14 x 12	16 x 12				
14" D	10 x 14	12 x 14	14 x 14	16 x 14	18 x 14			
16" D	10 x 16	12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	24 x 16
18" D		12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
20" D			14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
22" D				16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
24" D					18 x 24	20 x 24	22 x 24	24 x 24

*Odd sizes available in 1-inch increments.

Embrace part numbers

Cushion style	Stretch-Air	Comfort-Tek
Zero Elevation	EM-SZ- WWDD	EM-FZ- WWDD
Anti-Thrust	EM-SA- WWDD	EM-FA- WWDD

	Stretch-Air + GlideWear	Comfort-Tek + GlideWear
Zero Elevation	EM-SQZ- WWDD	EM-FQZ- WWDD
Anti-Thrust	EM-SQA- WWDD	EM-FQA- WWDD

Embrace ATI part numbers

Cushion style	Stretch-Air	Comfort-Tek
Zero Elevation	EM-SZM- WWDD	EM-FZM- WWDD
Anti-Thrust	EM-SAM- WWDD	EM-FAM- WWDD

	Stretch-Air + GlideWear	Comfort-Tek + GlideWear
Zero Elevation	EM-SQZM- WWDD	EM-FQZM- WWDD
Anti-Thrust	EM-SQAM- WWDD	EM-FQAM- WWDD

WWDD to be replaced | **WW** = width & **DD** = depth in inches

Weight capacity

14" wide & under	200 lbs
16" wide	250 lbs
18" wide	300 lbs
20" wide	350 lbs
22" wide & wider	400 lbs

M2 & M2 ATI



Features

- Promote postural alignment with tapered adductors to support the greater trochanters and deep lateral and medial contours to assist with leg positioning.
- Quadra gel pack and coccyx relief cutout can assist with pressure reduction on the pelvis and spine.
- Non-skid bottom fabric with available hook and loop attachments (HOOK-KIT-4).
- ATI base can be remolded to meet changes in both physical and functional needs.
- Warranty: 6 months cover / 5 years cushion

Dimensions 16" x 14" and smaller	Zero Elevation	Anti- Thrust
Seat thickness	3"	2.75"
Inferior thigh support thickness	2.5"	3"
Medial thigh support height	.75"	.75"
Lateral thigh support height	1"	1"
Pelvic contour height	1"	1"

Dimensions 16" x 16" and larger	Zero Elevation	Anti- Thrust
Seat thickness	4.25"	3.75"
Inferior thigh support thickness	3"	3.75"
Medial thigh support height	1.25"	1.25"
Lateral thigh Support height	1.5"	1.5"
Pelvic contour height	1.5"	1.25"

E2607 / E2608
Skin protection & positioning
E2624 / E2625
Adjustable skin protection & positioning

Available styles

Zero Elevation



Anti-Thrust



Standard sizes

	10" W	12" W	14" W	16" W	18" W	20" W	22" W	24" W
10" D	10 x 10	12 x 10	14 x 10					
12" D	10 x 12	12 x 12	14 x 12	16 x 12				
14" D	10 x 14	12 x 14	14 x 14	16 x 14	18 x 14			
16" D	10 x 16	12 x 16	14 x 16	16 x 16	18 x 16	20 x 16	22 x 16	24 x 16
18" D		12 x 18	14 x 18	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
20" D			14 x 20	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20
22" D				16 x 22	18 x 22	20 x 22	22 x 22	24 x 22
24" D					18 x 24	20 x 24	22 x 24	24 x 24

*Odd sizes available in 1-inch increments.

M2 part numbers

Cushion style	Stretch-Air	Comfort-Tek
Zero Elevation	M2-S- WWDD	M2-F- WWDD
Anti-Thrust	M2-SA- WWDD	M2-FA- WWDD

	Stretch-Air + GlideWear	Comfort-Tek + GlideWear
Zero Elevation	M2-SQ- WWDD	M2-FQ- WWDD
Anti-Thrust	M2-SQA- WWDD	M2-FQA- WWDD

M2 ATI part numbers

Cushion style	Stretch-Air	Comfort-Tek
Zero Elevation	M2-SM- WWDD	M2-FM- WWDD
Anti-Thrust	M2-SAM- WWDD	M2-FAM- WWDD

	Stretch-Air + GlideWear	Comfort-Tek + GlideWear
Zero Elevation	M2-SQM- WWDD	M2-FQM- WWDD
Anti-Thrust	M2-SQAM- WWDD	M2-FQAM- WWDD

WWDD to be replaced | **WW** = width & **DD** = depth in inches

Weight capacity

14" wide & under	200 lbs
16" wide	250 lbs
18" wide	300 lbs
20" wide	350 lbs
22" wide & wider	400 lbs

Comfort Custom Inception

E2609
Custom fabricated seat cushion



Dimensions

Width	10" -28" (1" increments)
Depth	10" - 28" (1" increments)

*Additional sizes may be available, contact Customer Experience for details.

Custom inception part number

C-C-**XW****XD**

X & **X** to be replaced | **XW** = width & **XD** = depth in inches

Tru-Shape®

E2609
Custom fabricated seat cushion



Features

- Custom molded seating solution with the width and depth designed to match an individuals shape and measurements.
- Evenly distribute weight to reduce pressure points and enhance skin protection.
- Can accommodate for asymmetries and unique positioning needs.
- No weight limit if appropriately sized to the individual.
- Warranty: 6 months cover/ 5 years cushion

ROHO Custom



Dimensions	
Cell heights	2.25" - 4.25"

E2609
Custom fabricated seat cushion

Features

- Customize the size, cells and layout
- No weight limit if appropriately sized to the individual.
- Warranty: 6 months cover / 3 years cushion

Custom part numbers

QSCUSTC	ROHO Quadtro Select (Custom cushion)
CUSTC	ROHO Single Compartment (Custom cushion)
ENHCUSTC	ROHO Enhancer (Custom cushion)
CSCUSTC	ROHO Contour Select (Custom cushion)

ROHO LTV Seat



Dimensions	
Cushion size	18" x 16"

Features

- A non-medical device that is lightweight and easy to fold making it perfect to take anywhere.
- Made of interconnected air cells that allow air to transfer from chambers for pressure redistribution.
- Available in one size with a quilted fabric cover option.
- Warranty: Cover 6 months / LTV cushion 1 year

LTV part number

LTVC

ROHO Pack-It



Dimensions	
Height	1.25"
Width	16.25"
Depth	9.5"
Average weight*	1.5 lb.

*Based on mid-sized cushioning device

Features

- Adjustable, air-filled, cellular-design.
- Neoprene surface intended for short-term cushioning needs.
- No weight limit yet the cushion must be properly sized to the individual.
- Warranty: 2 years

Pack-It part number

PACKIT

ROHO Mini-Max®



Dimensions	
Height	1.25"
Average weight*	2.25 lb.

*Based on mid-sized cushioning device

Standard sizes

	15" W	16" W	18" W	20" W
15" D	15 x 15			
16" D		16 x 16	18 x 16	
18" D			18 x 18	20 x 18
20" D				20 x 20

Features

- Adjustable, air-filled design intended as a supplemental seating surface.
- No weight limit if appropriately sized to the individual.
- Lightweight and easy to carry.
- Versatile use and excellent for sports, bathing, and travel.
- Made of non-permeable neoprene (not made with natural rubber latex).
- Warranty: 2 years

Mini-Max part number

MINI**WWDDC**

WWDD to be replaced | **WW** = width & **DD** = depth in inches

ROHO cushion sizing guide

ROHO neoprene standard series cushion sizing chart

Number of cells	Fits wheelchair size	Actual cushion size
5	10" (22cm - 26.5cm)	9.75"
6	11 - 12" (26.5cm - 31cm)	11.5"
7	13" (31cm - 35cm)	13.25"
8	14 - 15" (35cm - 39.5cm)	14.75"
9	16 - 17" (39.5cm - 44cm)	16.5"
10	18 - 19" (44cm - 48.5cm)	18.25"
11	20 - 21" (48.5cm - 52.5cm)	20"
12	22" (52.5cm - 57cm)	21.75"
13	23 - 24" (57cm - 61cm)	23.5"
14	25" (61cm - 65.5cm)	25"
15	26 - 27" (65.5cm - 70cm)	26.75"

ROHO neoprene select series (with ISOFLO control) cushion sizing chart

Number of cells	Fits wheelchair size	Actual cushion width**	Actual cushion depth**
5	10" (22cm - 26.5cm)	11"	10"
6	11 - 12" (26.5cm - 31cm)	12.75"	11.75"
7	13" (31cm - 35cm)	14.5"	13.5"
8	14 - 15" (35cm - 39.5cm)	16"	15.25"
9	16 - 17" (39.5cm - 44cm)	17.75"	17"
10	18 - 19" (44cm - 48.5cm)	19.5"	18.5"
11	20 - 21" (48.5cm - 52.5cm)	21.25"	20.25"
12	22" (52.5cm - 57cm)	23"	22"
13	23 - 24" (57cm - 61cm)	24.5"	23.75"
14	25" (61cm - 65.5cm)	26.25"	25"

**Measured from edge of cushion base material not edge of cells



Supports & accessories

Often overlooked, proper support accessories are a critical component to an effective seating & positioning system. Permobil specializes in providing unique, innovative accessories to ensure effective foot, arm, & amputee support.

BodiLink Head Support

Piano Taper (PT) hardware



MAX 8 shown

Gear Tooth (GT) hardware



Lateral shown

Lateral

Ideal for people who require lateral offset.

Mid-line

Ideal for people who require mid-line adjustment.

E1033
Removable hardware

MID 6

Mid-line adjustment offering 6" of reach

LAT 6

Lateral offset and mid-line adjustment offering 6" of reach.

MID 8

Mid-line adjustment offering 8" of reach.

LAT 8

Lateral offset and mid-line adjustment offering 8" of reach.

MAX 8

Maximum amount of adjustment offering 8" of reach.

Features

- GT hardware includes sintered-metal gear teeth joints with 15" increments, resulting in a strong positive-locking hardware PT hardware styles range in functionality and reach.
- Tool-less wing lock option for PT hardware for on-the-spot adjustments.
- Extra links available for creating a custom piece of hardware for user’s needs.
- 3" and 7" stem lengths available.

PT hardware part numbers

Size and shape	Part number
MID6	BL-HS-PT1RMMID6
LAT6	BL-HS-PT1RMLAT6
MID8	BL-HS-PT1RMMID8
LAT8	BL-HS-PT1RMLAT8
MAX8	BL-HS-PT1RMMAX8

PT flip down hardware part numbers

Size and shape	Part number
MID6	BL-HS-PT1FDMID6
LAT6	BL-HS-PT1FDLAT6
MID8	BL-HS-PT1FDMID8
LAT8	BL-HS-PT1FDLAT8
MAX8	BL-HS-PT1FDMAX8


GT hardware part numbers

Size and shape	Part number
Lateral	BL-HS-GT3RMLAT
Mid-line	BL-HS-GT3RMMID

Mounting hardware options



Posterior A or B mounting block

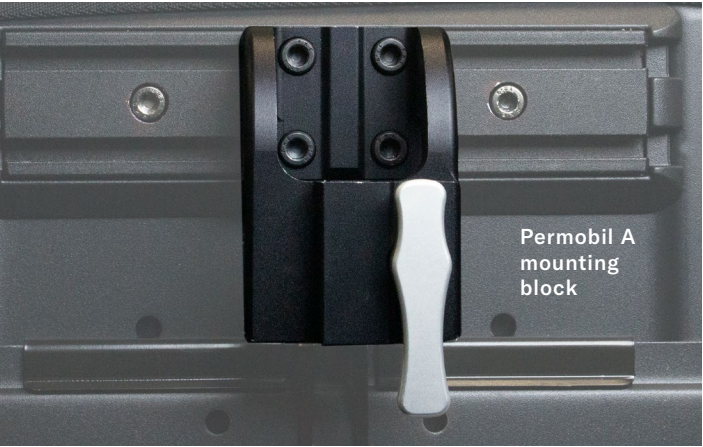

Permobil A mounting block



Relief bracket with mounting block (18W Shown)

Posterior A mounting block [\(HS-MB-POSTA\)](#)

For Embrace, Visco-Back, and Relief with Power Mount Only.

Posterior B mounting block [\(HS-MB-POSTB\)](#)

For OEM wood/t-nutback supports.

Anterior A mounting block [\(HS-MB-ANTA\)](#)

For AGILITY and other OEM back supports.

Quantum A mounting block [\(HS-MB-PW1A\)](#)

For Quantum Tru-Shape comfort back supports.

Permobil A mounting block [\(HS-MB-PW2A\)](#)

For Permobil Corpus 3G back supports (23"-28"H).

Permobil B mounting block [\(HS-MB-PW2B\)](#)

For Permobil Corpus 3G back supports (20"H).

Permobil C mounting block [\(HS-MB-PW2C\)](#)

For use with Permobil Direct Mount and Unitrak rail.

Relief bracket with mounting block

For Relief non-power mount only, includes mounting bracket and POSTA Block.

- 14W bracket with mounting block [\(HS-MB-AR14W\)](#)
- 16W bracket with mounting block [\(HS-MB-AR16W\)](#)
- 18W bracket with mounting block [\(HS-MB-AR18W\)](#)
- 20W bracket with mounting block [\(HS-MB-AR20W\)](#)

Premium pad supports

E0955
Cushioned head support pad



12W x 5L Soft Shown with GlideWear

Premium pad part numbers

Size	Pad firmness	Part number
6W x 3L	Soft	BL-HSP1-6W3L
6W x 3L	Extra soft	BL-HSP2-6W3L
8W x 4L	Soft	BL-HSP1-8W4L
8W x 4L	Extra soft	BL-HSP2-8W4L
10W x 5L	Soft	BL-HSP1-10W5L
10W x 5L	Extra soft	BL-HSP2-10W5L
12W x 5L	Soft	BL-HSP1-12W5L
12W x 5L	Extra soft	BL-HSP2-12W5L
14W x 5L	Soft	BL-HSP1-14W5L
14W x 5L	Extra soft	BL-HSP2-14W5L
18W x 5L	Soft	BL-HSP1-18W5L
18W x 5L	Extra soft	BL-HSP2-18W5L


Features

- Hinged wings allow for lateral contouring and support on 10W, 12W, 14W, and 18W pads.
- Soft and extra soft firmness foam available.
- 3 Cover Choices: Comfort-Tek, Stretch-Air, or GlideWear

Cover option

- Comfort-Tek ([COMFORT-TEK](#))
- Stretch-Air ([STRETCH-AIR](#))
- GlideWear ([GLIDEWEAR](#))

Basic pad supports

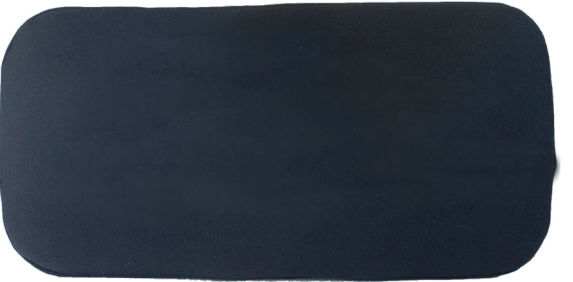

12W x 6L Shown with Stretch-Air

Basic pad part number

Size	Part number
6W x 5L	BL-HSP3-6W5L
8W x 5L	BL-HSP3-8W5L
10W x 6L	BL-HSP3-10W6L
12W x 6L	BL-HSP3-12W6L
14W x 6L	BL-HSP3-14W6L
18W x 6L	BL-HSP3-18W6L

Features

- Can be formed with moldable steel insert inside the pad.
- 2 Cover Choices: Comfort-Tek or Stretch-Air

Cover option

- Comfort-Tek ([COMFORT-TEK](#))
- Stretch-Air ([STRETCH-AIR](#))

BodiLink Lateral Trunk Support

E1034
Swing-away hardware

Gear Teeth 2 (GT2) series hardware



Swing-away
slot mount

Features

- 6 pad sizes offered in each premium & basic, moldable & rigid, supportive & active, giving 48 options to fit any user
- Spring-loaded joints for easy set-up and positioning.
- Unique wrist joint allows for pad rotation to match user shape.
- Additional links available for full customization.


Swing-away hardware part numbers

Style	Part number left	Part number right
Slot mounting	BL-LTS-GT2SAS-LH	BL-LTS-GT2SAS-RH
Relief slot mount	BL-LTS-GT2SAK-LH	BL-LTS-GT2SAK-RH
Embrace slot mount w/ winter/summer	BL-LTS-GT2SAWSAE-LH	BL-LTS-GT2SAWSAE-RH
Cane mounting	BL-LTS-GT2SAC-LH	BL-LTS-GT2SAC-RH
Quantum® power	BL-LTS-GT2SAPW1-LH	BL-LTS-GT2SAPW1-RH
Permobil® power	BL-LTS-GT2SAPW2-LH	BL-LTS-GT2SAPW2-RH
Track mount for Visco and Unite back	BL-LTS-GT2SAT-LH	BL-LTS-GT2SAT-RH

Fixed hardware part numbers

Style	Part number left	Part number right
Slot mounting	BL-LTS-GT2FXS-LH	BL-LTS-GT2FXS-RH
Relief slot mount	BL-LTS-GT2FXX-LH	BL-LTS-GT2FXX-RH
Embrace slot mount w/ winter/summer	BL-LTS-GT2FXWSAE-LH	BL-LTS-GT2FXWSAE-RH
Cane mounting	BL-LTS-GT2FXC-LH	BL-LTS-GT2FXC-RH
Quantum power	BL-LTS-GT2FXPW1-LH	BL-LTS-GT2FXPW1-RH
Permobil power	BL-LTS-GT2FXPW2-LH	BL-LTS-GT2FXPW2-RH
Track mount for Visco and Unite back	BL-LTS-GT2SFXT-LH	BL-LTS-GT2FXT-RH

Premium pad series

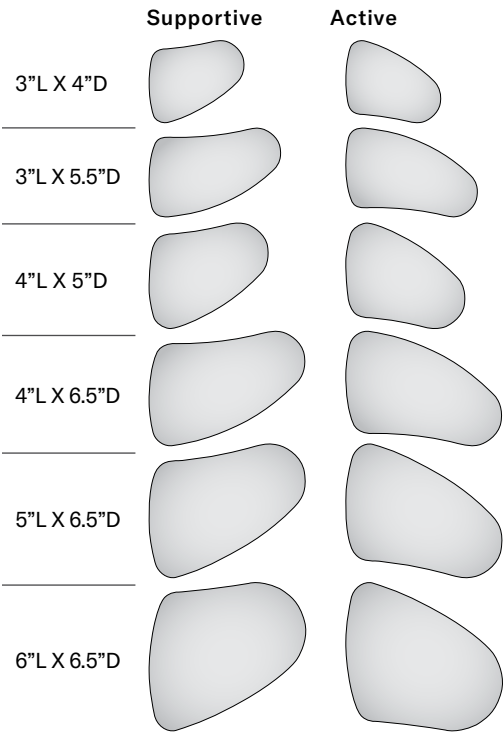

Features

- 7 Cover choices. Each available with added foam for an even softer feel.
- Anatomically designed pads in moldable or rigid style offer more natural positioning and greater surface contact.

Premium pad supportive part numbers

Size	Part number left	Part number right
3Lx4D	BL-LTSP2S-3L4D-LH	BL-LTSP2S-3L4D-RH
3Lx5.5D	BL-LTSP2S-3L5D-LH	BL-LTSP2S-3L5D-RH
4Lx5D	BL-LTSP2S-4L5D-LH	BL-LTSP2S-4L5D-RH
4Lx6.5D	BL-LTSP2S-4L6D-LH	BL-LTSP2S-4L6D-RH
5Lx6.5D	BL-LTSP2S-5L6D-LH	BL-LTSP2S-5L6D-RH
6Lx6.5D	BL-LTSP2S-6L6D-LH	BL-LTSP2S-6L6D-RH


Premium pad active part numbers

Size	Part number left	Part number right
3Lx4D	BL-LTSP2A-3L4D-LH	BL-LTSP2A-3L4D-RH
3Lx5.5D	BL-LTSP2A-3L5D-LH	BL-LTSP2A-3L5D-RH
4Lx5D	BL-LTSP2A-4L5D-LH	BL-LTSP2A-4L5D-RH
4Lx6.5D	BL-LTSP2A-4L6D-LH	BL-LTSP2A-4L6D-RH
5Lx6.5D	BL-LTSP2A-5L6D-LH	BL-LTSP2A-5L6D-RH
6Lx6.5D	BL-LTSP2A-6L6D-LH	BL-LTSP2A-6L6D-RH

Premium cover options

- No Cover ([LTS-P2-NO-CV](#))
 - Comfort-Tek ([LTS-P2-COMFORT](#))
 - Stretch-Air ([LTS-P2-STRETCH](#))
 - GlideWear ([LTS-P2-GLIDEWEAR](#))
- Comfort-Tek w/ Foam ([LTS-P2-COMFORT-FM](#))
 - Stretch-Air w/ Foam ([LTS-P2-STRETCH-FM](#))
 - GlideWear w/ Foam ([LTS-P2-GLIDEWEAR-FM](#))

Basic pad series

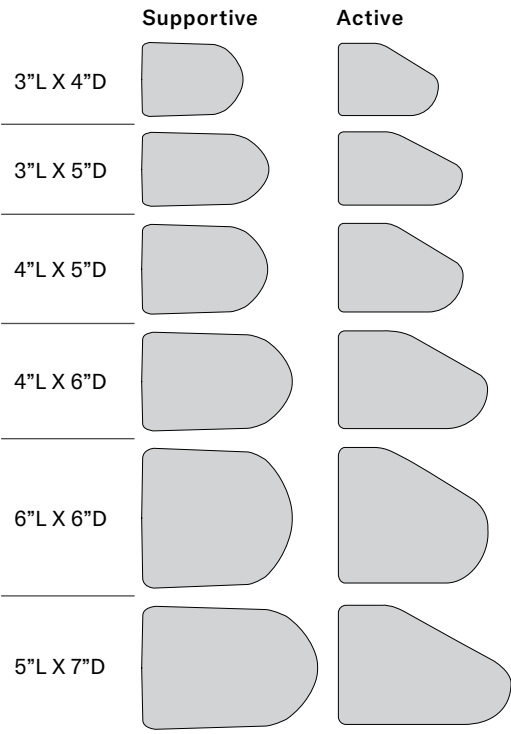

Features

- Pads available in moldable or rigid style.
- 2 Cover Choices: Comfort-Tek or Stretch-Air

Basic pad supportive part numbers

Size	Part number left	Part number right
3Lx4D	BL-LTSP1S-3L4D-LH	BL-LTSP1S-3L4D-RH
3Lx5D	BL-LTSP1S-3L5D-LH	BL-LTSP1S-3L5D-RH
4Lx5D	BL-LTSP1S-4L5D-LH	BL-LTSP1S-4L5D-RH
4Lx6D	BL-LTSP1S-4L6D-LH	BL-LTSP1S-4L6D-RH
5Lx6D	BL-LTSP1S-5L7D-LH	BL-LTSP1S-5L7D-RH
6Lx6D	BL-LTSP1S-6L6D-LH	BL-LTSP1S-6L6D-RH

Basic pad active part numbers

Size	Part number Left	Part number right
3Lx4D	BL-LTSP1A-3L4D-LH	BL-LTSP1A-3L4D-RH
3Lx5D	BL-LTSP1A-3L5D-LH	BL-LTSP1A-3L5D-RH
4Lx5D	BL-LTSP1A-4L5D-LH	BL-LTSP1A-4L5D-RH
4Lx6D	BL-LTSP1A-4L6D-LH	BL-LTSP1A-4L6D-RH
5Lx6D	BL-LTSP1A-5L7D-LH	BL-LTSP1A-5L7D-RH
6Lx6D	BL-LTSP1A-6L6D-LH	BL-LTSP1A-6L6D-RH

Basic cover options

- Comfort-Tek ([LTS-P1-COMFORT](#))
- Stretch-Air ([LTS-P1-STRETCH](#))

E0956
Lateral trunk support

BodiLink Lateral Pelvic/Thigh Support

E1034
Removable hardware

Telescoping Taper (TT) series



TT slot mount quick-release

Gear Tooth (GT) series



GT slot mount quick-release

TT slot mount removable hardware part numbers

Style	Part number left	Part number right
Small	BL-LPTS-TT1RMSL1-LH	BL-LPTS-TT1RMSL1-RH
Medium	BL-LPTS-TT1RMSL2-LH	BL-LPTS-TT1RMSL2-RH
Large	BL-LPTS-TT1RMSL3-LH	BL-LPTS-TT1RMSL3-RH

TT slot mount fixed hardware part numbers

Style	Part number left	Part number right
Small	BL-LPTS-TT1FXSL1-LH	BL-LPTS-TT1FXSL1-RH
Medium	BL-LPTS-TT1FXSL2-LH	BL-LPTS-TT1FXSL2-RH
Large	BL-LPTS-TT1FXSL3-LH	BL-LPTS-TT1FXSL3-RH

TT power mount removable hardware part numbers

Style	Part number left	Part number right
Small	BL-LPTS-TT1RMPWL1-LH	BL-LPTS-TT1RMPWL1-RH
Medium	BL-LPTS-TT1RMPWL2-LH	BL-LPTS-TT1RMPWL2-RH
Large	BL-LPTS-TT1RMPWL3-LH	BL-LPTS-TT1RMPWL3-RH

TT power mount fixed hardware part numbers

Style	Part number left	Part number right
Small	BL-LPTS-TT1FXPWL1-LH	BL-LPTS-TT1FXPWL1-RH
Medium	BL-LPTS-TT1FXPWL2-LH	BL-LPTS-TT1FXPWL2-RH
Large	BL-LPTS-TT1FXPWL3-LH	BL-LPTS-TT1FXPWL3-RH

GT slot mount removable hardware part numbers

Part number left	Part number right
BL-LPTS-GT2RMS-LH	BL-LPTS-GT2RMS-RH


GT slot mount fixed hardware part numbers

Part number left	Part number right
BL-LPTS-GT2FXS-LH	BL-LPTS-GT2FXS-RH

Features

- TT Power Mount Hardware mounts to: Permobil, Quantum, Quickie, Rovi, and Avid power seating systems.
- TT Hardware has precision telescoping height and angle adjustment.
- Extra links available for GT and long extension arm available for TT hardware for customization.
- Fixed or removable options with high visibility release levers. TT featuring a one-handed release lever.

TT hardware size vs. cushion thickness

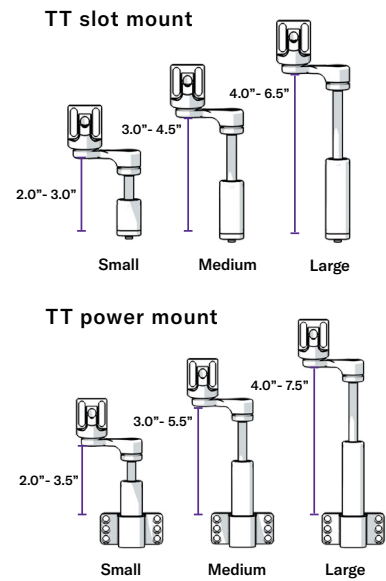

Features

- 4 premium pad sizes and 10 basic pad sizes.
- Basic pads available Zero-Elevation
- Premium pads are constructed of molded foam to reduce breakdown.
- 3 cover choices: Comfort-Tek, Stretch-Air, and GlideWear

Premium pad series

Premium pads	
3.5”L X 4”D	
3.5”L X 8”D	
3.5”L X 12”D	
5.5”L X 6”D	

Premium pad part numbers

Size	Part number left	Part number right
3.5Lx4D	BL-LPTSP2Z-3L4D-LH	BL-LPTSP2Z-3L4D-RH
3.5Lx8D	BL-LPTSP2Z-3L8D-LH	BL-LPTSP2Z-3L8D-RH
3.5Lx12D	BL-LPTSP2Z-3L12D-LH	BL-LPTSP2Z-3L12D-RH
5.5Lx6D	BL-LPTSP2Z-5L6D-LH	BL-LPTSP2Z-5L6D-RH

Premium cover options


- Comfort-Tek ([LPTS-P2-COMFORT-FM](#))
- Stretch-Air ([LPTS-P2-STRETCH-FM](#))
- GlideWear ([LPTS-P2-GLIDE-FM](#))

Basic cover options

- Comfort-Tek ([LPTS-P1-COMFORT](#))
- Stretch-Air ([LPTS-P1-STRETCH](#))
- GlideWear ([LPTS-P1-GLIDE](#))

BodiLink Medial Knee/Thigh Supports are also available, please visit [permobil.com](#) for information.

Basic pad series

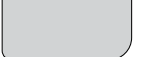Basic Zero-Elevation pads	
3”L X 4”D	
3”L X 5”D	
4”L X 4”D	
4”L X 6”D	
5”L X 7”D	
4”L X 8”D	
4”L X 10”D	
4”L X 12”D	
4”L X 14”D	
4”L X 16”D	

Basic pad Zero-Elevation part numbers

Size	Part number left	Part number right
3Lx4D	BL-LPTSP1Z-3L4D-LH	BL-LPTSP1Z-3L4D-RH
3Lx5D	BL-LPTSP1Z-3L5D-LH	BL-LPTSP1Z-3L5D-RH
4Lx4D	BL-LPTSP1Z-4L4D-LH	BL-LPTSP1Z-4L4D-RH
4Lx6D	BL-LPTSP1Z-4L6D-LH	BL-LPTSP1Z-4L6D-RH
4Lx8D	BL-LPTSP1Z-4L8D-LH	BL-LPTSP1Z-4L8D-RH
4Lx10D	BL-LPTSP1Z-4L10D-LH	BL-LPTSP1Z-4L10D-RH
4Lx12D	BL-LPTSP1Z-4L12D-LH	BL-LPTSP1Z-4L12D-RH
4Lx14D	BL-LPTSP1Z-4L14D-LH	BL-LPTSP1Z-4L14D-RH
4Lx16D	BL-LPTSP1Z-4L16D-LH	BL-LPTSP1Z-4L16D-RH
5Lx7D	BL-LPTSP1Z-5L7D-LH	BL-LPTSP1Z-5L7D-RH

Comfort Foot Single



IST - Inferior support thickness
L - Posterior support length

Features

- Unique heel cut-out eliminates contact with heel, prevents pressure points, and reduces risk of skin breakdown.
- Optional gel pack in foot support base reduces pressure, shear, and heat build-up.
- Medium density, high resiliency foam construction provides maximum comfort.
- Non-skid fabric on foot support helps maintain foot position.
- Tapered bottom accommodates new style footplates.
- Fits on both left and right foot.
- Standard with Comfort-Tek cover.

Comfort Foot Single

Style	W x D x L	Part number	Gel part number
Standard	7 x 10 x 8	FTS	FTSG
Tall	7 x 10 x 12	FTST	FTSTG
Bariatric	10 x 10 x 8	FTSB	FTSBG

Inferior support pad (base only)

Style	W x D	Part number	Gel part number
Standard	7 x 10	FTS-BTM	FTSG-BTM
Bariatric	10 x 10	FTSB-BTM	FTSBG-BTM

Comfort Foot Single elevating

Style	W x D x L	Part number
Standard	7 x 10 x 8	FTSE
Tall	7 x 10 x 12	FTSTE

Hardware

Clamp-on connection: Easy on, easy off connection gives users a versatile attachment.

Bolt-on connection: Bolt-on connection works with all foot plates.



Style	Inferior support thickness (IST)
Standard	2”
W/gel	2.5”

Style	Posterior support thickness
Standard	1.75”
W/gel	1.75”

Posterior support pad (top portion only)

Style	W x D x L	Part number
Standard top portion only	7 x 10 x 8	FTS-TOP
Tall top portion only	7 x 10 x 12	FTST-TOP
Bariatric top portion only	10 x 10 x 8	FTSB-TOP
Elevating top portion only	7 x 10 x 8	FTSE-TOP
Tall elevating top portion only	7 x 10 x 12	FTSTE-TOP

Reimbursement for items billed using K0108 is not guaranteed, as coverage is at the carrier’s discretion.

Comfort Foot Double



Standard



Tall



E0954
Footbox

Features

- Unique heel cut-out eliminates contact with heel, prevents pressure points, and reduces risk of skin breakdown.
- Optional medial lower limb separator keeps feet and legs positioned properly.
- Non-skid fabric on foot support helps maintain foot position.
- Durable securing buckle and bolt-on option comes standard with all models.
- Optional gel pack in foot support reduces pressure.
- Comes standard with Comfort-Tek cover.

Posterior support pad (top portion only)

WC size	W x D x L	W/separator	W/out separator
16"	15.5" x 10.5" x 8"	FTD-16-TOP	FTDNS-16-TOP
18"	17.5" x 10.5" x 8"	FTD-18-TOP	FTDNS-18-TOP
20"	19.5" x 10.5" x 8"	FTD-20-TOP	FTDNS-20-TOP
22"	21.5" x 10.5" x 8"	FTD-22-TOP	FTDNS-22-TOP
24"	23.5" x 10.5" x 8"	FTD-24-TOP	FTDNS-24-TOP

Tall posterior support pad (top portion only)

WC size	W x D x L	W/separator	W/out separator
16"	15.5" x 10.5" x 12.5"	FTDT-16-TOP	FTDTNS-16-TOP
18"	17.5" x 10.5" x 12.5"	FTDT-18-TOP	FTDTNS-18-TOP
20"	19.5" x 10.5" x 12.5"	FTDT-20-TOP	FTDTNS-20-TOP
22"	21.5" x 10.5" x 12.5"	FTDT-22-TOP	FTDTNS-22-TOP
24"	23.5" x 10.5" x 12.5"	FTDT-24-TOP	FTDTNS-24-TOP

Inferior support pad (base only)

WC size	W x D	Part number	Gel part number
16"	15.5" x 10.5"	FTD-16-BTM	FTDG-16-BTM
18"	17.5" x 10.5"	FTD-18-BTM	FTDG-18-BTM
20"	19.5" x 10.5"	FTD-20-BTM	FTDG-20-BTM
22"	21.5" x 10.5"	FTD-22-BTM	FTDG-22-BTM
24"	23.5" x 10.5"	FTD-24-BTM	FTDG-24-BTM

Bolt-on installation

Bolt-on installation with securing buckle allows Comfort Foot Double to swing away for transfers. Product can be installed with securing buckles only.

Style	Inferior support thickness (IST)
Standard	2"
W/gel	2.5"

Comfort Foot Double w/leg separator

WC size	W x D x L	Part number	Gel part number
16"	15.5" x 10.5" x 8"	FTD-16	FTDG-16
18"	17.5" x 10.5" x 8"	FTD-18	FTDG-18
20"	19.5" x 10.5" x 8"	FTD-20	FTDG-20
22"	21.5" x 10.5" x 8"	FTD-22	FTDG-22
24"	23.5" x 10.5" x 8"	FTD-24	FTDG-24

Comfort Foot Double w/out leg separator

WC size	W x D x L	Part number	Gel part number
16"	15.5" x 10.5" x 8"	FTDNS-16	FTDNSG-16
18"	17.5" x 10.5" x 8"	FTDNS-18	FTDNSG-18
20"	19.5" x 10.5" x 8"	FTDNS-20	FTDNSG-20
22"	21.5" x 10.5" x 8"	FTDNS-22	FTDNSG-22
24"	23.5" x 10.5" x 8"	FTDNS-24	FTDNSG-24

Comfort Foot Double Tall w/leg separator

WC size	W x D x L	Part number	Gel part number
16"	15.5" x 10.5" x 12.5"	FTDT-16	FTDTG-16
18"	17.5" x 10.5" x 12.5"	FTDT-18	FTDTG-18
20"	19.5" x 10.5" x 12.5"	FTDT-20	FTDTG-20
22"	21.5" x 10.5" x 12.5"	FTDT-22	FTDTG-22
24"	23.5" x 10.5" x 12.5"	FTDT-24	FTDTG-24

Comfort Foot Double Tall w/out leg separator


WC size	W x D x L	Part number	Gel part number
16"	15.5" x 10.5" x 12.5"	FTDTNS-16	FTDTNSG-16
18"	17.5" x 10.5" x 12.5"	FTDTNS-18	FTDTNSG-18
20"	19.5" x 10.5" x 12.5"	FTDTNS-20	FTDTNSG-20
22"	21.5" x 10.5" x 12.5"	FTDTNS-22	FTDTNSG-22
24"	23.5" x 10.5" x 12.5"	FTDTNS-24	FTDTNSG-24

Flexion Footbox

Features

- Single Footbox fits on both left and right side and is available in four sizes: 5", 6", 7", & 8" wide.
- Double Footbox fits on four wheelchair sizes: 14", 16", 18", and 20" wide.
- Easy angle adjustment from 90° to 160° with a single screw.
- Durable aluminum shell.
- After drilling mounting holes in footbox, bolt on to chair for permanent placement.
- All surfaces and edges are foam padded for added comfort.
- Comes standard with Comfort-Tek cover.

E0954
Footbox

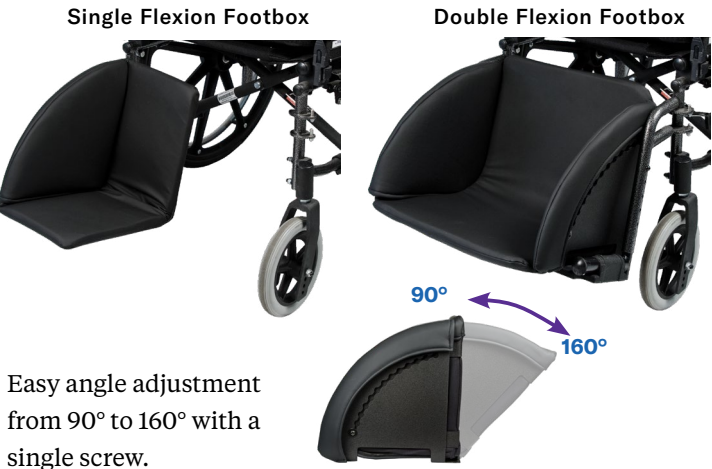

Easy angle adjustment from 90° to 160° with a single screw.

Single Footbox

W x D x L	Part number
5.25" X 11" X 11"	FBS-5
6.25" X 11" X 11"	FBS-6
7.25" X 11" X 11"	FBS-7
8.25" X 11" X 11"	FBS-8

Double Footbox

W x D x L	WC size	Part number
12" X 11" X 11"	14"	FBD-14
14" X 11" X 11"	16"	FBD-16
16" X 11" X 11"	18"	FBD-18
18" X 11" X 11"	20"	FBD-20

Calf protector



E0995
Calf protector

Features

- Provides lower extremity support.
- Buckle attachment makes transfers easy.
- Available in both standard and rigid options.
- Standard with Comfort-Tek cover.

Calf protector

WC size	W x L	Thickness	Part number	Rigid part number
16"	14" X 8"	0.75"	70A	71A
18"	16" X 8"	0.75"	70B	71B
20"	18" X 8"	0.75"	70C	71C
22"	20" X 8"	0.75"	70D	71D

Reimbursement for items billed using K0108 is not guaranteed, as coverage is at the carrier's discretion.

Quick-Release Amp

E1020
Residual limb support system



Features

- Residual limb support may be used with a solid seat pan.
- Installs onto either left or right side of solid seat pan.
- Toggle lever allows positioning at multiple angles.
- Pad has adduction to help center residual limb.
- Foam construction for comfort yet a firmer base layer for support.
- Comes standard with Comfort-Tek cover.
- TB3 parts designed for use on Quantum Rehab Tru-Balance® 3 chairs.

10" Quick-Release Amp part numbers

W x D	Thickness	Extension kit
8" X 10"	3"	AMPKIT810
9" X 10"	3"	AMPKIT910
10" X 10"	3"	AMPKIT1010
11" X 10"	3"	AMPKIT1110
12" X 10"	3"	AMPKIT1210

10" TB3 Quick-Release Amp part numbers

W x D	Thickness	Extension kit
8" X 10"	3"	AMPTQR810
9" X 10"	3"	AMPTQR910
10" X 10"	3"	AMPTQR1010
11" X 10"	3"	AMPTQR1110
12" X 10"	3"	AMPTQR1210

Quick-Release Amp hardware only

Part number
AMPQHARD

Quick-Release Amp hardware

Fast and simple removal with the quick release assembly.



Angle adjustment

Toggle lever allows positioning at multiple angles.



14" Quick-Release Amp part numbers

W x D	Thickness	Extension kit
8" X 14"	3"	AMPKIT814
9" X 14"	3"	AMPKIT914
10" X 14"	3"	AMPKIT1014
11" X 14"	3"	AMPKIT1114
12" X 14"	3"	AMPKIT1214

14" TB3 Quick-Release Amp part numbers

W x D	Thickness	Extension kit
8" X 14"	3"	AMPTQR814
9" X 14"	3"	AMPTQR914
10" X 14"	3"	AMPTQR1014
11" X 14"	3"	AMPTQR1114
12" X 14"	3"	AMPTQR1214

Quick-Release Amp hardware only

Part number
AMPSTQRHARD

Swing-Away Amp

E1020
Residual limb support system




Swing-Away Amp hardware

Hardware allows for fast and simple swing-away and removal.

Pad can be easily removed from and reinstalled onto wheelchair.

Hardware comes in aluminum, aluminum large, or steel.

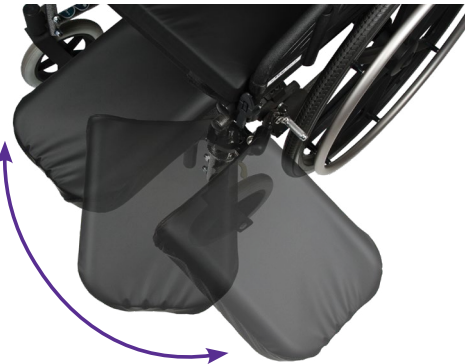

Angle adjustment

Toggle lever allows positioning at multiple angles.

Up to 4" of height adjustment to accommodate varying seat cushion thicknesses.

Pad can be adjusted 3" in depth to accommodate wheelchair depth conversion kits.



Features

- Provides a fully adjustable residual limb support.
- Allows individuals to transfer in and out of the wheelchair easily.
- Aluminum and Steel fits tubing sizes from 7/8" to 1".
- Aluminum large fits tubing sizes from 1" to 1 1/4".
- Pad has adduction to help center residual limb.
- Inferior support thickness is 3".
- Comes standard with Comfort-Tek cover.

Swing-Away Amp aluminum

Width	10" depth	14" depth
8"W	AMPSA810	AMPSA814
9"W	AMPSA910	AMPSA914
10"W	AMPSA1010	AMPSA1014
11"W	AMPSA1110	AMPSA1114
12"W	AMPSA1210	AMPSA1214

Swing-Away Amp aluminum large

Width	10" depth	14" depth
8"W	AMPSALG810	AMPSALG814
9"W	AMPSALG910	AMPSALG914
10"W	AMPSALG1010	AMPSALG1014
11"W	AMPSALG1110	AMPSALG1114
12"W	AMPSALG1210	AMPSALG1214

Swing-Away Amp steel

Width	10" depth	14" depth
8"W	AMPSAS810	AMPSAS814
9"W	AMPSAS910	AMPSAS914
10"W	AMPSAS1010	AMPSAS1014
11"W	AMPSAS1110	AMPSAS1114
12"W	AMPSAS1210	AMPSAS1214

Swing-Away Amp hardware only

Type	Part number
Aluminum	AMPSAHARD
Aluminum large	AMPSALGHARD
Steel	AMPSASHARD

Comfort Arm

E2209
Arm support
E1028
Arm support hardware



Standard Arm



Low Arm



Tall Rail Arm




Molded Arm

Features

- Comfort Arm design provides finger and hand support.
- Comfortable straps for securing can be adjusted or removed.
- Provides support for contractures and tone.
- Moldable option can be used to elevate fingers and hand.
- Standard, Low and Tall Rail arm come standard with Comfort Tek cover.

Comfort Arm dimensions

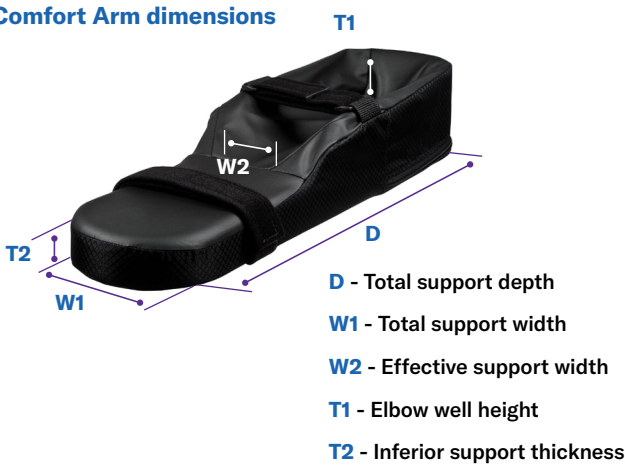

Elevating articulating depth-adjustable (EAD) hardware

Allows arm support to be elevated and articulated with a quick release lever.

Arm support can be adjusted for a depth of up to 2.5” and a width of up to 4”.



Elevating



Articulating



Depth-adjustable

Interlock bracket

The interlock bracket locks down to 3/4” up to 1” tubing with a two-bolt connection.

Standard with fixed or swing-away hardware option.

Swing-away hardware

The rotating base system allows the Comfort Arm to be rotated side-to-side and adjusted front to back 7”.

Slide on hardware

Slides arm support over stock wheelchair arm pad.

TB3 hardware

Arm support with track mount hardware for TB3 style chairs.

Comfort Arm supports

Standard Arm supports

D x W1	W2	T1	T2	Side	Part number
22” X 6”	4.0”	2.75”	1.5”	Right	2002R-PAD
22” x 6”	4.0”	2.75”	1.5”	Left	2002L-PAD
14” x 6”	4.0”	2.75”	1.5”	Right	2002R-14-PAD
14” x 6”	4.0”	2.75”	1.5”	Left	2002L-14-PAD
17” x 6”	4.0”	2.75”	1.5”	Right	2002R-17-PAD
17” x 6”	4.0”	2.75”	1.5”	Left	2002L-17-PAD

Low Arm supports

D x W1	W2	T1	T2	Side	Part number
22” x 6”	3.0”	n/a	1.25”	Right	2004R-PAD
22” x 6”	3.0”	n/a	1.25”	Left	2004L-PAD
14” x 6”	3.0”	n/a	1.25”	Right	2004R-14-PAD
14” x 6”	3.0”	n/a	1.25”	Left	2004L-14-PAD
17” x 6”	3.0”	n/a	1.25”	Right	2004R-17-PAD
17” x 6”	3.0”	n/a	1.25”	Left	2004L-17-PAD

Molded Arm supports

Depth	Side	Part number
22”	Right	2300R-PAD
22”	Left	2300L-PAD
14”	Right	2300R-14-PAD
14”	Left	2300L-14-PAD

Standard Arm moldable supports

D x W1	W2	T1	T2	Side	Part number
22” x 6”	4.0”	2.75”	1.5”	Right	2002R-M-PAD
22” x 6”	4.0”	2.75”	1.5”	Left	2002L-M-PAD

Low Arm moldable supports

D x W1	W2	T1	T2	Side	Part number
22” x 6”	3.0”	n/a	1.25”	Right	2004R-M-PAD
22” x 6”	3.0”	n/a	1.25”	Left	2004L-M-PAD

Molded Arm for EAD hardware

Depth	Side	Part number
22”	Right	2300PR-PAD
22”	Left	2300PL-PAD
14”	Right	2300PR-14-PAD
14”	Left	2300PL-14-PAD

Tall Rail Arm supports

D x W1	W2	T1	T2	Side	Part number
22” x 7”	4.0”	2.75”	1.5”	Right	2007R-PAD
22” x 7”	4.0”	2.75”	1.5”	Left	2007L-PAD

Tall Rail Arm moldable supports

D x W1	W2	T1	T2	Side	Part number
22” x 7”	4.0”	2.75”	1.5”	Right	2007R-M-PAD
22” x 7”	4.0”	2.75”	1.5”	Left	2007L-M-PAD

Comfort Arm hardware

Hardware type	Side	Part number
Fixed interlock	Right	1000R-HARD
Fixed interlock	Left	1000L-HARD
Swing-away	Right	2000R-HARD
Swing -away	Left	2000L-HARD
EAD	Right	3000R-HARD
EAD	Left	3000L-HARD
Slide on	Right	2100R-HARD
Slide on	Left	2100L-HARD
TB3 V3 track bottom mount	Left	ARMBTM-P3L-HARD
TB3 V3 track bottom mount	Right	ARMBTM-P3R-HARD

Durasoft Full Tray

E0950
Tray



Dimensions

Thickness	1"
Width	24"
Depth	18"

Features

- Rigid insert is covered by high resiliency foam and easily cleaned Comfort-tek cover.
- Quick attachment to full length or desk length arms.
- Stability due to webbing or slide mount attachment.
- Forward adjustment for larger users.
- Easily moved for transfers.

Durasoft Full Tray part number

Mounting	Part number
Webbing	75W
Side mount	75K-SM

Flip-up Half Tray

E0950
Tray



Dimensions

Thickness	1"
Width	12"
Depth	21"

Features

- Rigid insert is covered by high resiliency foam and easily cleaned Comfort-Tek cover.
- Interlock attachment version provides quick attachment to full length or desk length arms.
- Slide-on version provides quick-on, quick-off installation.
- Easily flipped 270 degrees out of the way for transfers.

Flip-up Half Tray part number

Mounting	Side	Part number
Interlock	Right	761FR
Interlock	Left	761FL
Slide on	Right	761SFR
Slide on	Left	761SFL

Durasoft Half Tray

E0950
Tray



Dimensions

Thickness	1.5" - 3"
Width	12"
Depth	21"

Features

- Raised outside edge keeps arm positioned correctly.
- Easily swiveled out of the way for transfers or table use.
- Interlock attachment system provides quick attachment to full length or desk length arms.

Durasoft Half Tray part number

Option	Side	Part number
Tray and hardware	Right	762R
Tray and hardware	Left	762L
Tray only	Right	762R-TRAY
Tray only	Left	762L-TRAY
Hardware only	Right	700R-HARD
Hardware only	Left	700L-HARD

Rigid Insert

E0992
Solid seat support base



Features

- Decreases the “hammock” effect in outstretched wheelchair upholstery and provides extra support for bariatric users.
- Designed to fit all Comfort Company cushions.
- Made of ABS plastic.
- Comes in 3/16” thickness.

Rigid Insert part number

RIDWWDD
WWDD to be replaced WW = width & DD = depth in inches
Sizes Available: 10"W x 10"D - 28"W x 22"D

Removable Solid Seat Pan

E2231
Solid seat support base



Features

- Includes a slotted aluminum pan and attaching hardware for 7/8” or 1” tubing.
- Comes with two different cross bars to accommodate various wheelchair frame types.
- Lateral thigh supports can be easily attached.

Removable Solid Seat Pan part number
RSSP-1-WWDD
WWDD to be replaced | WW = width & DD = depth in inches
Sizes Available: 10”W x 10”D - 28”W x 22”D

Incontinent Proof Liner

E2619
Replacement cover



Features

- Provides extra incontinence protection for seating surfaces.
- The liner may be used as an extra cover when the outer cover is being washed.
- Cushion style and size must be specified at time of order.


Please contact CE for part numbers.

Replacement cover

E2619
Replacement cover




GlideWear™
TECHNOLOGY
Available with
GlideWear technology.

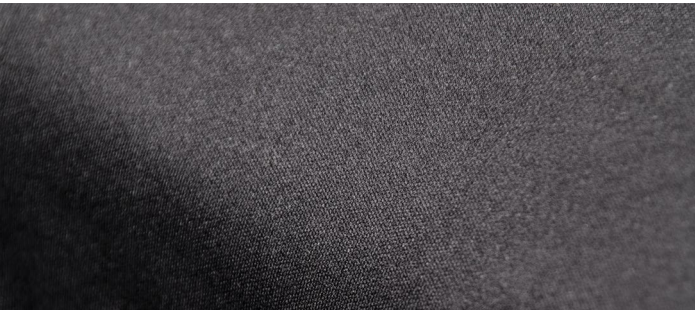

Stretch-Air

Stretch-Air is designed to be airy and breathable with multi-directional stretch for optimal immersion into the cushion.

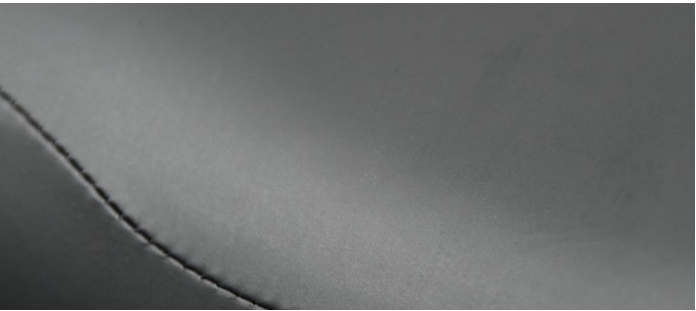

Comfort-Tek

Comfort-Tek is smooth, strong and soft with multi-directional stretch for immersion into the cushion. Designed for infection control and as an incontinence barrier, it is easily cleaned with common disinfectants and is latex-free.

Replacement covers are available for seating surfaces and back supports.

ROHO Planar Solid Seat Insert

E0992Manual wheelchair accessory, solid seat insert



Standard stock sizes								
	12" W	14" W	15" W	16" W	17" W	18" W	20" W	24" W
12" D	12 x 12							
14" D	12 x 14	14 x 14						
15" D		14 x 15	15 x 15					
16" D		14 x 16	15 x 16	16 x 16				
17" D			15 x 17	16 x 17	17 x 17			
18" D				16 x 18	17 x 18	18 x 18		
20" D				16 x 20		18 x 20	20 x 20	
22" D						18 x 22	20 x 22	
24" D						18 x 24	20 x 24	24 x 24

Features

- Provides a solid base for use with most ROHO cushions.
- Lightweight and easily modified.
- Prevent the hammocking effect of wheelchair sling upholstery.
- Easily fits under the cushion, inside the ROHO cushion cover.
- Made of polyurethane coated birch plywood, the 3/8"/ (1 cm).

Planar Solid Seat Insert part number
SIWWDD
WWDD to be replaced WW = width & DD = depth in inches

Maximum size is 24" x 24"

Comfort Privacy Shield




Features

- Connects with hook and loop attachment to back canes for coverage between a back support and cushion.

Privacy Shield part number	
BS-PS-1015W	For chairs 10"W - 15"W
BS-PS-1621W	For chairs 16"W - 21"W
BS-PS-2225W	For chairs 22"W - 25"W
BS-PS-2628W	For chairs 26"W - 28"W

ROHO Adaptor Pad

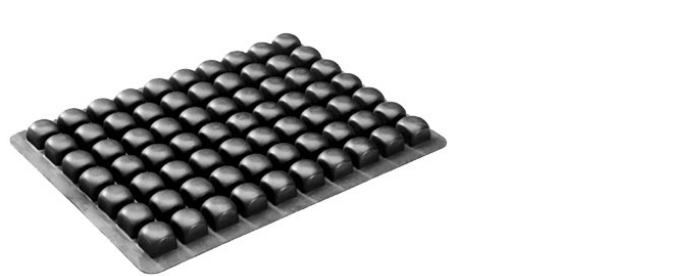

Adaptor Pad part number	
ADAPTOR	
Dimensions	
Cell height	1"
Width	9.25"
Depth	13"

Features

- The self-contained air cells can be used as an interface on virtually any flat or contoured surface to protect skin and soft tissues.
- Can be cut to fit various shapes such as sides of wheelchairs, trays, head and footrests, inside of helmets, or for any surfaces that may need protection.
- Easily cleaned and disinfected.
- No weight limit if appropriately sized to the individual.

ROHO Contour Base

E0992Manual wheelchair accessory, solid seat insert



Standard sizes			
	13" W	15" W	17" W
13" D	13 x 13		
15" D		15 x 15	17 x 15
17" D			17 x 17

Features

- Molded base that provides a solid foundation for use with ROHO cushions.
- Easily fits under the cushion inside the ROHO cushion cover.
- Lightweight.
- Reduces the effects of sling upholstery.
- Assists in posture and positioning.

Contour Base part numbers
BASE1313
BASE1515
BASE1715
BASE1717

ROHO Recliner System



Recliner System part number

RECLINERCUSTC

Features

- Adjustable, air-filled, cellular-design, support surface that utilizes Dry Floatation® technology and is intended to conform to an individual's seated shape to provide localized skin/soft tissue protection.
- Custom recliner cushion and back system.
- Made to fit per order.
- There is no weight limit, yet the system must be properly sized to the individual.
- Warranty: 2 years

ROHO Shower/Commode Seat



Dimensions

Cell height	2.25"
Width	14.75" - 18.25"
Depth	14.75" - 18.25"

Standard sizes

	14.75" W	16.5" W	18.25" W
14.75" D	14.75 x 14.75		
16.5" D	14.75 x 16.5	16.5 x 16.5	
18.25" D		16.5 x 18.25	18.25 x 18.25

Features

- Designed to fit shower or commode chairs.
- Available in 5 standard sizes; also available in custom sizes.
- Does not include cover.
- Easy to clean and disinfect.
- No weight limit if sized appropriate to individual.
- Warranty: 2 years

Shower/Commode Seat part numbers

SEATCMD88LP	SEATCMD910LP
SEATCMD89LP	SEATCMD1010LP
SEATCMD99LP	SEATCMDCUST

ROHO Toilet Seat cushion



Dimensions

Width	7.75"
Depth	15.75"

Features

- An adjustable, air-filled cellular-design, neoprene surface intended for short-term cushioning needs.
- Easy to attach and remove – straps and snaps securely fasten the cushion to the toilet seat.
- Multi-height cells: outer cells 2.25” tall, inner cells 1” tall.
- Does not include cover.
- Clean with soap and water.
- No weight limit.
- Warranty: 2 years

Toilet Seat cushion part number

SEATTOILET

References

Weight conversion

Pounds	Kilograms	Pounds	Kilograms
100 lbs	45.4 kg	350 lbs	158.8 kg
150 lbs	68.0 kg	400 lbs	181.4 kg
200 lbs	90.7 kg	450 lbs	204.1 kg
250 lbs	113.4 kg	550 lbs	249.5 kg
300 lbs	136.1 kg	750 lbs	340.2 kg
330 lbs	149.7 kg		

Size conversion

Inches	Centimeters	Inches	Centimeters
10"	25.4 cm	22"	55.9 cm
12"	30.5 cm	24"	61.0 cm
14"	35.6 cm	26"	66.0 cm
16"	40.6 cm	28"	71.1cm
18"	45.7 cm	30"	76.2 cm
20"	50.8 cm		

© 2025 Permobil, Inc. Products may be covered by one or more U.S. and foreign patents and trademarks, including but not limited to: THE ADAPTOR®, AGILITY®, AIRLITE®, CONTOUR SELECT®, DRY FLOATATION®, HIGH PROFILE®, Hybrid Elite SR®, LOW PROFILE®, MINI-MAX®, QUADTRO SELECT®, ROHO®, Sensor Ready™, Smart Check®, Comfort-Tek™, Compass®, Curve®, QUADROGEL®, and Elements®.

Product information is changed as needed; current product information is available at permobil.com.

ROHO, Inc. and The Comfort Companies, LLC are part of Permobil, Inc.

GlideWear™ is a trademark of Tamarack Habilitation Technologies, Inc. BOA® is a trademark of BOA Technology Inc. Zytel® is a trademark of Dupont Polymers, Inc. Vicair®, SmartCells® and Vector® are trademarks of Vicair B.V.



Permobil offers mobility, seating and positioning options no matter what your age, size, shape or ability. We understand that every individual is unique in their needs. That's why we are proud to offer a Solution for every Body. Scan the QR code to learn more.